




Fisherbrand™ Products - The Ultimate Collection V3

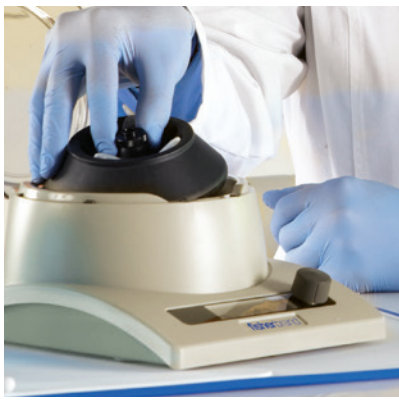

Contents

| | |
|-------------------------------------------------------------------------------------------------|----|
| Workspace Essentials <i>Centrifuges, Stirrers, Lab Mats, Notebooks & Drying Racks</i> | 4 |
| Liquid Handling <i>Pipettes, Reservoirs & Chillers, Pipette Racks & Stands</i> | 20 |
| Racks <i>Tubes, Microtubes, Freezers & Storage Racks</i> | 27 |
| Histology and Microscopy <i>Slide Boxes & Mailers, Staining Trays</i> | 47 |
| Sample and Blood Collection <i>Counters, Collection Trays & Storage, Transport Boxes</i> | 52 |
| Cryogenics <i>Boxes, Vial Canes & Sleeves, Storage Rack Systems</i> | 57 |
| Weighing <i>Papers, Dishes, Boats, Scoops & Spatulas</i> | 61 |
| Microbiology <i>Lead Rings, Cell Spreaders & Inoculating Picks</i> | 66 |
| Organisation <i>Bins, Dispensers, Containers & Organisers</i> | 67 |
| Sealing and Labelling <i>Parafilm & Foil Dispensers, Heat Sealers, Label Stations</i> | 73 |

Featured New Products



Clinical Centrifuge
Pg. 13



Triple Instrument Bundle
Pg. 15



PCR Bundle
Pg. 16



32 Place PCR Rack
Pg. 45



Colony Counter Bundle
Pg. 53



Cryogenic & Vaccine Vial Grippers
Pg. 58

Fisherbrand™ Centrifuge/Stirrer

Multi-functional Lab Solution

Saves Valuable Benchtop Space

Perfect for Small Research and Education Labs

Holds
5 mL
Tubes



High-Speed Mini Centrifuge



Two
Units
in One

5
Year
Warranty



Magnetic Stirrer

EASY SWITCH OVER
NO
TOOLS REQUIRED

High-speed Centrifuge is easy to use and designed to deliver maximum speed (12,500 rpm/9,783 x g) and flexibility with rapid acceleration and fast controlled breaking.

Ideal for most protocols requiring fast spins including PCR, micro- filter cell separation and HPLC. Versatility to handle 0.2, 0.5, 1.5, 2.0 and 5.0 mL tube sizes.

Intuitive One-Dial Programming

- Set speed and time with ease:
- Adjustable speed from 500 to 12,500 rpm
- Speed setting for RPM or RCF
- Adjustable time from 30 seconds to 60 minutes

Safe to Use

Safety features help ensure it will not operate when the lid is open, a rotor is jammed, or the unit is out of balance.

High-performance Magnetic Stirrer is designed using high quality rare earth magnets creating a robust magnetic coupling that significantly reduces the chance of spin-outs.

Perfect for work with sensitive samples and protocols up to 3 L requiring aggressive mixing for full dilution.

Easy to Use

Set stir mode, speed and time with easy one-dial programming:

- Select Mode: clockwise, counter-clockwise and oscillate
- Adjustable speed from 50 to 2,500 rpm
- Continuous and adjustable time settings from 30 seconds to 60 minutes

Durable Accessories

- Includes reusable, easy-to-clean silicone lid mats to stabilize vessels and prevent spills and two stir bars (PTFE-coated rare-earth magnets, 40x12 mm)

| Cat. No. | Colour | L x W x H cm | Qty |
|----------|--------|------------------|-----|
| 15909350 | Gray | 23.5 x 18.5 x 12 | 1 |

| Replacement parts | | | |
|-------------------|-----------------------------|-------------------------------------|---|
| 16380886 | High Speed Centrifuge Rotor | 12x1.5/2.0 mL Tube | |
| 16390886 | High Speed Centrifuge Rotor | 6 x 5 mL tube | |
| 16300896 | Stirrer Rotor | | 1 |
| 16320896 | Silicone Lid Mat | | |
| 16310896 | Stir Bar | 40 x 12 mm for Centrifuge / Stirrer | |
| 16370886 | Tube Adapters | 0.2 mL, 0.5 mL | 6 |

Patent Information

Community Design No. 003515337-0001-0002
Euro. Pat. App. 3 246 088
US Design Patent No. D814045
Utility Patent Pending



| Package Includes | | |
|--------------------|---------------------------------|-------------------------|
| Centrifuge/Stirrer | Stirrer Rotor | 12 Tube Adapters 0.2 mL |
| 5 mL Tube Rotor | 2 Stir Bars | 12 Tube Adapters 0.5 mL |
| 1.5 mL Tube Rotor | 2 Silicone Mats for Stirrer Lid | |

1 Low Voltage, Double-Insulated Power Adapter with 4 Interchangeable Plugs



| Specifications: Centrifuge | |
|----------------------------|----------------------------------------------------------------------------------------------------------------|
| Speed Range | 500 to 999 rpm (Increments of 1); 1,000 to 5,000 rpm (Increments of 5); 5,000 to 12,500 rpm (Increments of 25) |
| RCF | 16 to 9,783 x g (Max RCF Calculated) |
| Run Time Ranges | 30 Seconds to 60 Minutes (5 Second Increments) |

| Specifications: Magnetic Stirrer | |
|----------------------------------|----------------------------------------------------------------------------|
| Speed Range | 50 to 999 rpm (increments of 1); 1,000 to 2,500 rpm (increments of 5) |
| Run Time Ranges | Stirrer Mode: 30 Seconds to 60 Minutes or Continuous (5 Second Increments) |
| | Oscillating Stirrer Mode: 1 Second to 60 Minutes or Continuous |
| | Oscillating Period: 30 Seconds to 15 Minutes |
| Capacity | 3 Liters Plus. Diameter of Vessel not to Exceed 130 mm |

Fisherbrand™ High-Speed Mini-Centrifuge

Maximum speed with minimal space

This compact centrifuge is ideal for protocols calling for fast spins with a higher RCF requirement. Made of durable ABS with UV-Resistant coating.



User Friendly

Intuitive push-button programming and bright LCD screen make it easy to use.

Quiet Operation

Brush-less motor enables quiet operation and low vibrations.

Safety Features

Lid will not open while rotor is turning and clear polycarbonate allows user to see rotor movement for added visual safety. Also includes automatic shut-off feature for sample imbalance and tilt.



| Cat. No. | Description | L x W x H cm | Qty |
|----------|---------------------------------------------------|--------------------|-----|
| 12972041 | Mini Centrifuge 100-240 VAC, 50/60 Hz | 20.3 x 17.1 x 11.4 | |
| 12992041 | Universal Power Cord Set | | |
| 12902051 | Rotor 12 x 1.5/2.0 mL Tube (with Cover and Knob) | | |
| 12912051 | Rotor 4 x 8 strip 0.2 mL Tube (with Knob) | | 1 |
| 12922051 | Rotor Knob | | |
| 12932051 | Rotor Cover (for 12 x 1.5/2.0 mL Tube Rotor Only) | | |
| 12982041 | Power Adapter | | |

Package includes:

- High-Speed Mini Centrifuge
- 12-Place 1.5/2.0 mL Tube Rotor
- 32-Place 0.2 mL Tube Rotor (singles, or 4 x 8 Strip)
- 12 Tube Adapters 0.2 mL
- 12 Tube Adapters 0.5 mL
- Universal Double-Insulated Power Adapter (110–240 V 50/60 Hz) with 4 Interchangeable Plugs



Tool-free rotor changes

Specifications:

Speed: 500 to 12,500 rpm
 Timer Range: 1 to 30 minutes
 RCF: 15 to 9,783 x g
 Five-year warranty



Fisherbrand™ Micro & Conical Centrifuge Tubes

Pure black for light sensitive samples

Microcentrifuge Tubes: 1.5 mL (a)

- Attached flat caps for easy one-handed use
- Frosted writing area
- Syringe pierceable cap

Conical Centrifuge Tubes: 15 and 50 mL (b, c)

- Sterile with lot-to-lot traceability
- Writing area and flat polyethylene lids for easy labelling



Certified DNase, RNase and Pyrogen free

Autoclavable

Microcentrifuge Tubes

| | Cat. No. | Colour | Vol. mL | Outer Dimensions mm | Closed Height mm | Graduations mL | Max RCF Rating | Packaging | Qty |
|---|----------|--------|---------|---------------------|------------------|------------------|----------------|-----------|-----|
| a | 15386548 | Black | 1.5 | 13.4 | 40.6 | 0.1, 0.5, 1, 1.5 | 18,000 x g | Bulk | 500 |

Conical Centrifuge Tubes

| | Cat. No. | Colour | Vol. mL | Outer Dimensions mm | Closed Height mm | Graduations mL | Max RCF Rating | Packaging | Qty |
|---|----------|--------|---------|---------------------|------------------|----------------|----------------|-------------------|-----|
| b | 15396548 | Black | 15 | 17 | 120 | 0.5 | 9,400 x g | Bulk Wrap Sterile | 500 |
| c | 15306558 | Black | 50 | 29 | 114 | 5.0 | 9,400 x g | Bulk Wrap Sterile | 500 |

Fisherbrand™ Mini-Centrifuge

Mini but Mighty - ideal for small runs

The Mini Centrifuge is ideal for microfilter cell separations and HPLC samples, with adapters to accommodate all major microtube and PCR tube styles (1.5/2.0, 0.5, 0.2 mL) for customised use.



Easy to Use

Tool-free SnapSpin rotors and quick-spin options make the centrifuges fast and easy to use.

Ergonomic Design

- The lid fits in the palm of the hand for one-handed locking, unlocking, and quick tube removal. The lid height accommodates taller lidded and skirted tubes
- The ABS plastic housing and polycarbonate plastic lid make cleaning easy

Safety Measures

- Fitted with a braking function, the unit quickly stops when powered off or when lid is opened
- A weighted rubberised base eliminates potential movement and dampens vibration and noise



Accessories and Storage Case Included



Specifications:

Maximum Speed: Approx. 6,000 rpm (fully loaded)
 Maximum RCF: 2,000 x g
 Power: 100–240 VAC 50/60 Hz 0.6A, Approx. 7.3 W
 Five-year warranty



Package includes:

- Mini-Centrifuge
- 6-Place 1.5/2.0 mL Tube Rotor
- 16-Place 0.2 mL Tube Rotor (singles, or 2 x 8 strip)
- 6 Tube Adapters 0.2 mL
- 6 Tube Adapters 0.5 mL
- 4 O-Rings
- Storage Case for Rotors and Adapters
- Tube Rack
- Low Voltage, Double-Insulated Power Adapter with 4 Interchangeable Plugs



SnapSpin Rotors

change out with the push of a button. No tools are needed

| Cat. No. | Colour | L x W x H cm | Qty |
|----------|--------|--------------------|-----|
| 15358266 | Gray | 15.3 x 12.8 x 10.9 | 1 |

| Replacement Parts | Qty |
|---------------------------------------------|-----|
| 12883332 SnapSpin Rotor 6 x 1.5/2.0 mL Tube | 1 |
| 12893332 Snapspin Rotor 2 x 8 Strip 0.2 mL | 1 |

Fisherbrand™ Portable Centrifuge Kit

Spin samples on the go



Specially Designed Inserts

seat the centrifuge, all its accessories, and includes a tube rack



12 Volt Car Adapter
for on-the-go testing

The Portable Centrifuge Kit offers a convenient way to spin down and prepare samples outside the laboratory environment.

Immediate Field Processing

- Process samples in remote areas before deterioration can take place
- Allows for immediate on-site testing for faster results
- Made of durable ABS plastic for easy cleaning in field environments

Easy to Use

- Centrifuge is equipped with tool-less SnapSpin rotors for easy change out
- Simple operation with lock/unlock activation

Portable Convenience

- All the kit items fit in a Transport Box for portability
- The handled box has an air-tight seal for optimal sample protection



| Cat. No. | Description | Qty |
|----------|--------------------------------------------------------------------------------------|-----|
| 15817944 | Portable Centrifugation Kit, 100-240 VAC, 50/60 Hz Universal Plug + 12 V Car Adapter | 1 |
| 15827944 | 12V Mini Centrifuge Car Adapter* | 1 |

*Adapter fits Fisherbrand Mini Centrifuge units manufactured since Jan 2016

Package includes:

Mini Centrifuge

- 6-Place 1.5/2.0 mL SnapSpin Tube Rotor
- 16-Place 0.2 mL SnapSpin Tube Rotor (singles or 2 x 8 strip)
- 6 Tube Adapters (0.2 mL)
- 6 Tube Adapters (0.5 mL)
- 4 Lid O-Rings
- Low Voltage, Double-Insulated Power Adapter with 4 Interchangeable Plugs



Transport Box

- 16-Place Tube Rack Insert
- Transport Box Dividers
- Mini Centrifuge Cradle
- Portable 12 V Mini Centrifuge Car Adapter
- Five-year warranty



Fisherbrand™ Magnetic Stirrer

Strong, motor-less and maintenance free


Superior Functional Design

- Advanced magnetic coil technology drives the spinning of the stir bar maintaining continuous, constant speed for reproducible results
- Magnetic coupling is optimized by using alternating current (AC) as its driving force strong enough to stir viscous solutions
- Motor-less operation generates minimal heat ideal for stirring temperature sensitive samples and environments for long periods
- With no moving parts, produces less vibration and quieter operation over traditional magnetic stirrers

Intuitive Operation

- Multiple stirring modes are easy to program
- Bright digital display precisely controls time and speed
- Quickly meets stirring speeds with rapid acceleration and braking (less than 10 seconds)

Protective silicone anti-skid plate cover for improved safety

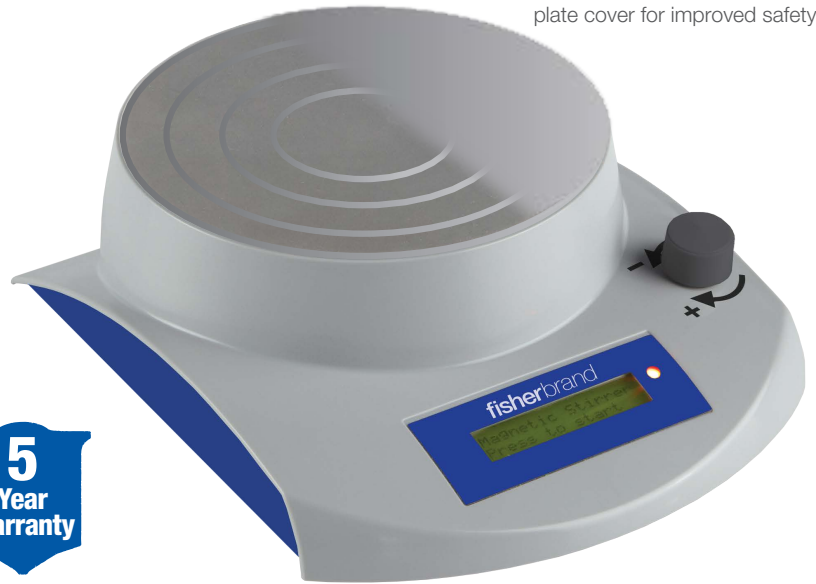

5
Year
Warranty

NEW

Durable ABS housing is corrosion resistant and easy to clean



Specifications

| | |
|----------------|-----------------------------------------------------------------------|
| Speed Range | 50–2,000 rpm variable, 1 rpm increments |
| Run Time Range | 30 seconds to 60 minutes, variable, 1 second increments or continuous |
| Stirrer Modes | Single clockwise |
| | Single counter clockwise |
| | Auto reverse rotation 1 to 99 cycles 50% duty |
| Capacity | 3 L plus |
| Mixing Surface | 165 mm diameter, silicone mat 160 mm diameter |



Package Includes

- Rare Earth, NdFeB Magnetic Induction Stirrer
- 2 Stir Bars 35 mm x 12 mm Cylindrical Shape
- 1 Silicone Anti-Skid Plate Cover
- 1 Low Voltage, Double Insulated Power Adapter (100-240 VAC, 50/60 Hz) with 4 Interchangeable Plugs



Magnetic Stirrer

| Cat. No. | Color | L x W x H | | Qty |
|----------|-----------|---------------|-------------------|------|
| | | in | cm | |
| 17891672 | Grey/Blue | 9.7 x 7.7 x 3 | 24.5 x 19.5 x 7.5 | 1 ea |



**5
Year
Warranty**



NEW

Fisherbrand™ Mini Stirrer

Affordable solution, ultimate function

The Mini Magnetic Stirrer is an affordable, dependable solution ideal for general research, academic, life sciences, industrial and clinical lab settings

Durable Design

- Constructed of chemical resistant ABS and PC blend
- Compact footprint saves valuable space on benchtop
- Lightweight for easy transport

Ultimate Function

- Powerful magnetic coupling ensures accurate alignment, reducing spin-outs during acceleration or deceleration
- Performs routine stirring applications with low vibration and quiet operation



Specifications

Speed Range 200–2,200 rpm (Infinitely adjustable)

Capacity Up to 3 L




Package Includes

Mini Magnetic Stirrer

2 Stir Bars 35 mm x 9 mm, Octagonal Shape

1 Low Voltage, Double Insulated Power Adapter (100-240 VAC, 50/60 Hz) with 4 Interchangeable Plugs

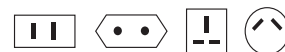

| Cat. No. | Color | L x W x H in | L x W x H cm | Qty |
|----------|---------------|------------------|-------------------|------|
| 17881672 | Grey/Blue ■ ■ | 7.2 x 5.2 x 1.85 | 18.3 x 13.2 x 4.7 | 1 ea |

U.S. and E.U. Utility Patent Pending
U.S. Design Patent No. D814,045



Fisherbrand™ Vortexer

Multipurpose, maximum convenience

NEW

Infinity Adjustable Speed Rotary Dial
Provides your choice of speed control



Innovative Design

- One mixing head does it all with this sleek, modern and compact designed vortex mixer
- Hold slot configurations for multiple sizes of tubes, vessels and assay plates
- Center pad for quick hand-held mixing

Easy to Use

- Control switch with three operation selections (intermittent, off, continuous)
- Variable Speed Rotary Dial provides pre-set speed control selection during shaking mode

Durable Construction

- ABS material resist most commonly used chemicals and is easy to clean
- Stable elastomer feet dampen vibration and prevent movement during operation



Center pad area for hand held mixing



Assay plate x 1 (to be used with retention cords)



50 mL conical tubes x 2 (horizontally mounted)



15 mL conical tubes x 2 (horizontally mounted)



5 mL tubes x 2 (horizontally mounted in the same retention space as 15 mL)



1.5/2 mL microfuge tubes x 4 (vertically mounted)



0.5 mL microfuge tubes x 6 (vertically mounted)



0.2 mL microfuge tubes (strip tubes) x 24 (vertically mounted)

| Cat. No. | Power Supply |
|----------|-----------------------------------------|
| 17801682 | 220V – 240V, 50-60 Hz, 44 Watts, CE, UK |



Specifications

| | |
|--------------------|------------------------------------------------------------------------------|
| Speed Range: | To 3,000 rpm, (shaking to vortexing) |
| Orbit: | 3.7 mm |
| Control: | 3 position slider — (top) intermittent, (middle) off, (bottom) continuous |
| Weight: | 8.66 Lb/3.93 Kg |
| Dimensions: | 171w x 190d x 197h mm (6.67w x 7.41d x 7.68h in.)(with mixing head in place) |
| Five-year warranty | |
| Quantity: | 1 ea |

Fisherbrand™ Mini Vortexer

small footprint with maximum power

The Mini Vortexer offers immediate vortex shaking for a variety of tube sizes. A simple touch start operation enables exceptional mixing performance with one hand

Robust Design

- Low profile design and small footprint saves valuable space
- Made of ABS and silicone materials that resist most commonly used chemicals
- Stable platform and elastomer feet ensure unit stays stable during operation

Powerful Function

- Rapid vortex action with low activation force
- Quiet operation and low vibration

5
Year
Warranty

NEW

Flexible

Can be used with tubes ranging from 2.0 mL to 50 mL

Non-slip Feet

ensures the unit stays in place while mixing



U.S. and E.U. Utility Patent Pending
U.S. Design Patent No. D814,045

| Item No. | Color | L x W x H in | L x W x H cm |
|----------|-----------|-----------------|-----------------|
| 17871672 | Grey/Blue | 5.1 x 4.3 x 3 | 13 x 11 x 7.5 |

Package Includes

Mini Vortexer

2 Silicone Caps

1 Low Voltage, Double Insulated Power Adapter (100-240 VAC, 50/60 Hz) with 4 Interchangeable Plugs



Specifications

| | |
|----------|------------------|
| Speed | 4,500 rpm |
| Orbit | 6 mm |
| Capacity | Up to 50 mL tube |
| Quantity | 1 ea, 8 ea/cs |





Fisherbrand™ Orbital Shaker

Shaker for a vast range of protocols

Versatile Design

- Create your own matrix for vessels of all shapes and sizes with Hi-Lo Platform and Bungee Cords
- Instantly adjusts to weight load changes and unbalanced loads with no-stop electronic motor feedback allowing for truly reproducible results
- Durable ABS construction suitable for use in incubators (excluding CO₂) and cold rooms within the stated temperature range

High-Featured Operation

- Suitable for wide range of protocols including very low start speed for sensitive sample agitation to high speeds for vigorous applications
- Soft start and shut down feature reduces spilling from full vessels
- Designed to handle heavy loads up to 10 lbs (4.5 kg)
- Overload protection shut-off prevents overloading



Digital Orbital Shaker

| Item No. | L x W x H | | Platform Interior W x D | | Qty |
|----------|-------------------|------------------|-------------------------|-------------|------|
| | in | cm | in | cm | |
| 17811682 | 12.4 x 7.8 x 12.2 | 31.4 x 19.7 x 31 | 11 x 11 | 27.9 x 27.9 | 1 ea |

Specifications

| | |
|----------------|-----------------------------------------------------------------|
| Speed Range | 20–300 rpm, 1 rpm increments |
| Orbit Diameter | 19 mm (0.8") |
| Run Time Range | 0 to 74 hours and 59 minutes; 1 second increments or continuous |
| Maximum Load | 4.5 kg (10 lbs) |



Package Includes

Digital Orbital Shaker

Skid-Free Rubber Platform Mat

8 Bungee Cord Set


Universal Hi-Lo Platform with Soft Rubber Edge Bumpers

1 Low Voltage, Double Insulated Power Adapter (100–240 VAC, 50/60 Hz) with 4 Interchangeable Plugs





Tube Adapter Sleeve

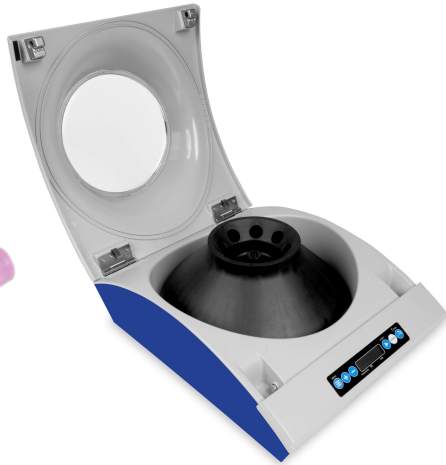

IVDR Designation

Class A, Rule 5a (Article 47, Annex VIII) per IVDF 2017/746

Applications include:

- Clinical chemistry
- Clinical microbiology
- Blood Collection Prep
- Blood Banking
- Immunology
- Pathology



Fisherbrand™ Clinical Centrifuge

Powerful clinical centrifugation in compact bench-top size

Superb Function

- Powerful, Adjustable Speed (500 to 6,000 RPM) can be set in RPM or RCF
- Programmable with memory for up to 8 programs
- Adjustable time (30 seconds to 60 minutes or continuous)
- Rapid acceleration and gentle electronic braking (deceleration time of <60 seconds)
- Digitally controlled brushless motor for accurate speed control and no motor maintenance required

Optimal Design

- Easy to read Digital back-lit display
- Tool-less rotor removal for convenient loading and cleaning
- Clear lid/Viewing porthole to allow non-contact speed verification by non-contact tachometer
- Bench-top Size - compact footprint for workspace use


Built-in Safety

- Safety features help ensure it will not operate when the lid is open, unit is jammed or unit is out of balance
- Manual lid unlock allows users to open lid in case of a power failure and prevents loss of samples



Rotor fits a multitude of tube sizes with the included tube adapters and spacers

- 3 - 10 mL Blood Collection Tubes
- 5 mL MicroCentrifuge Tubes
- 15 mL Conical Centrifuge Tubes

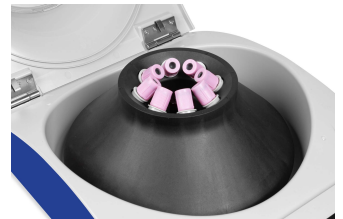

Package Includes

Clinical Centrifuge 8 Tube Adapters 8-Place Rotor 24 Tube Spacers

1 Low Voltage, Double-Insulated Power Adapter with 4 Interchangeable Plugs



Specifications: Clinical Centrifuge

| | |
|-----------------|--------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| Speed Range | 500 - 6,000 rpm / 33 x 4,668 x g (max RCFcalculated) with typical 15 mL Tube |
| Capacity | 8-place 15 mL tube rotor: 120 mL total maximum liquid volume |
| Run Time Ranges | 30 seconds - 60 minutes or continuous |
| Certifications | CE, UKCA, RoHS3, WEEE, IVDR* *Class A, Rule 5a (Article 47, Annex VIII) per IVDF 2017/746 |
| Weight | 4.5 kg (9.93 lb) |
| Electrical | Power Adapter Rating Input: 100-240 VAC, 2.5A, 50/60 Hz / Output: 12 VDC, 12.5 A, 150W Centrifuge Rating Input: 12 VDC / Current (Under Load): Approx. 7.7 A / Power: Approx. 93W |
| Operating Temp. | Range 2°C to 40°C |

| Cat. No. | Colour | L x W x H in | L x W x H cm | Qty |
|----------|-----------|-------------------|--------------------|-----|
| 17801727 | Grey/Blue | 14.4 x 11.5 x 5.8 | 36.6 x 29.3 x 14.8 | 1 |

Reusable Lab Mats

Sustainable, Safe Solutions

Reusable lab mats are an essential lab solution to help keep benchtops organised and clean. They help facilitate workflow processes and are ideal for every type of laboratory.

Save Money & Eliminate Waste

Since the mats are reusable, they reduce cost and waste making them a more environmentally friendly option than paper mat alternatives.

Adds Stability & Safety

These mats offer a large working surface that is anti-skid to keep fragile glassware and expensive equipment in place.

The mats surface is also noise dampening, reducing sound from the vibration of common benchtop equipment.

Made from durable food-grade silicone, they are chemically resistant and have a raised edge to help contain spills.

Easy to Clean & Store

Between uses, the mats can be easily wiped clean using any standard disinfectant or placed in a glassware washer.

Fisherbrand™ Workstation Lab Mat

Safely facilitate workflow processes

Gridded Surface

for tube staging, preventing movement of objects



Side One: combination gridded and smooth surface

Side Two: smooth surface only

- Protect countertops from hot items up to 200°C (392°F)
- Lip around the edge can hold approx. 250 mL of liquid
- Small footprint, 12 x 14.9 x 0.1 in. (30.5 x 37.9 x 0.2 cm)
- **Autoclavable** at 121° for 15 minutes

| Cat. No. | Colour | L x W x H cm | Qty |
|----------|--------|-------------------|-----|
| 16361648 | Blue | 30.5 x 37.9 x 0.2 | 1 |

Fisherbrand™ Silicone Lab Mat

Large working surface

The lab mat can be used with items such as centrifuges, vortex mixers, hot plates, stirrers, pipettes and so much more.

Mat features designated magnetic area to keep stir bars in place.

- Protect countertops from hot items up to 200°C (392°F)
- **Autoclavable** at 121° for 15 minutes

Large Working Surface

accommodates benchtop equipment, weigh boats, stir bars, scales, etc.

Raised Edges

retain up to 300 mL of spilled liquids



Reversible

Blue and Gray

| Cat. No. | Side #1 Colour | Side #2 Colour | L x W x H cm | Qty |
|----------|----------------|----------------|--------------|-----|
| 15837944 | Blue | Blue/Gray | 35 x 60 x .2 | 1 |

Fisherbrand™ Triple Instrument Bundle

Whether you need to spin-down, shake, or stir, this bundle has you covered.

Get the full functionality of three machines in two pieces of equipment! Complete with a Fisherbrand Silicone Lab Mat to protect your workspace from spills and wear.



Fisherbrand Centrifuge/Stirrer

The industry's first high-speed centrifuge (12,500 rpm) and magnetic stirrer (3 L), combined into one unit. Made durable and chemically resistant with easy cleaning and maintenance-free operation. Its compact and low-profile design optimizes limited bench space and minimizes clutter in the lab.

Fisherbrand Mini Vortexer

Small footprint Vortexer offering maximum vortex shaking for a variety of tube sizes. A simple touch start operation allows for immediate sample mixing with one hand and automatic stop when tube is lifted.

Fisherbrand Silicone Lab Mat

This FDA food-grade silicone labmat is suitable for various laboratory settings, including cold rooms, clean rooms, and sterile suites, where high sterility is essential. Effectively safeguards countertops from temperatures up to 200°C (392° F) and demonstrates general chemical resistance to moderate oxidizing chemicals.

- **Autoclavable** at 121° for 15 minutes



Reversible
Blue and Gray

| Cat. No. | Description | Qty |
|----------|---------------------------------------|-----|
| 17862966 | Fisherbrand Triple Instrument Bundle* | 1 |

*For detailed spec information, please refer to product pages: Centrifuge/Stirrer (4), Mini Vortexer (11), Silicone Lab Mat (13)

Fisherbrand™ PCR Bundle

Experience seamless PCR sampling, preparation, and storage with this specially curated bundle.



Fisherbrand Mini Centrifuge

This compact centrifuge offers user-friendly features such as tool-free SnapSpin rotors, quick-spin options, a braking function for quick stops, and a weighted rubberized base for stability and noise reduction.



SnapSpin Rotors
change out with the push of a button. No tools are needed!

Fisherbrand Mini Vortexer

Small footprint Vortexer offering maximum vortex shaking for a variety of tube sizes. A simple touch start operation allows for immediate sample mixing with one hand and automatic stop when tube is lifted.

Fisherbrand PCR Cooler

Ice-free cooling system for PCR preparation and storage. Cools industry-standard PCR tubes, strips, and plates with efficiency. Maintains the entire 96-Well plate at a temperature of 0°C for over an hour and keeps samples below 5°C for up to three hours. Vibrant color change signals to user when temperatures exceed 6°C.

Fisherbrand 96 Place Storage Rack Low Temp

These Low Temperature PCR Racks have a rotating sliding lid that doubles as a stand, providing convenience and stability. The translucent lid is specially designed for rapid cooling and quick sample identification. With ergonomic finger grips and a robust ABS base, these racks ensure secure handling and sample safety. Available in multiple colors for lab coding.

Fisherbrand Silicone Lab Mat

This FDA food-grade silicone labmat is suitable for various laboratory settings, including cold rooms, clean rooms, and sterile suites, where high sterility is essential. Effectively safeguards countertops from temperatures up to 200°C (392° F) and demonstrates general chemical resistance to moderate oxidizing chemicals.

- **Autoclavable** at 121° for 15 minutes

Reversible
Blue and Gray



Fisherbrand Mini Centrifuge (16617645)

| Colour | L x W x H in | L x W x H cm |
|--------|-----------------|--------------------|
| Grey | 6 x 5 x 4.3 | 15.3 x 12.8 x 10.9 |

Specifications

| | |
|--------------------|----------------------------------|
| Speed Range | Approx. 6,000 rpm (fully loaded) |
| RCF | 2,000 x g |



Package Includes

| | | |
|----------------------------------------------------------------------------------------------------|------------------------------------------------------|------------------------|
| Mini Centrifuge | 6-Place 1.5/2.0 mL Tube Rotor | 6 Tube Adapters 0.2 mL |
| 6 Tube Adapters 0.5 mL | 4 Lid O-Rings | Tube Rack |
| Storage Case for Rotors and Adapters | 16-Place 0.2 mL Tube Rotor (singles, or 2 x 8 strip) | |
| 1 Low Voltage, Double Insulated Power Adapter (100-240 VAC, 50/60 Hz) with 4 Interchangeable Plugs | | |



Fisherbrand Mini Vortexer (17871672)

| Colour | L x W x H in | L x W x H cm |
|-----------|-----------------|-----------------|
| Grey/Blue | 5.1 x 4.3 x 3 | 13 x 11 x 7.5 |

Specifications

| | |
|-----------------|------------------|
| Speed | 4,500 rpm |
| Orbit | 6 mm |
| Capacity | Up to 50 mL tube |



Package Includes

| | |
|----------------------------------------------------------------------------------------------------|-----------------|
| Mini Vortexer | 2 Silicone Caps |
| 1 Low Voltage, Double Insulated Power Adapter (100-240 VAC, 50/60 Hz) with 4 Interchangeable Plugs | |



Fisherbrand PCR Cooler (16320412)

| Warm Colour | Cold Colour | Wells | Tube Sizes mL | L x W x H in | L x W x H cm |
|-------------|-------------|-------|------------------|-----------------|-----------------|
| Pink | Purple | 96 | 0.2 Assay Plates | 4.3 x 5.8 x 1.9 | 11 x 14.8 x 4.8 |
| Light Blue | Dark Blue | | | | |



Fisherbrand Silicone Lab Mat (15837944)

| Colour | L x W x H in | L x W x H cm |
|-----------|-----------------|-----------------|
| Grey/Blue | 14 x 23.5 x 0.1 | 35 x 60 x 0.2 |



96 Place Storage Rack Low Temp (15758177)

| Colour | Tubes mL | Wells | Well Dia. mm | L x W x H in | L x W x H cm |
|-----------|-------------|-------|-----------------|-----------------|------------------|
| Assorted* | 0.2 | 96 | 5.2 | 4.9 x 3.6 x 1.5 | 12.5 x 9.2 x 3.7 |

*Assorted colours (Purple, Blue, Green, Orange, Yellow)






Keep Information Safe
with a durable hardbound,
water-resistant cover

Easy Identification
with a foil stamped title

Marker Ribbon
to keep your place

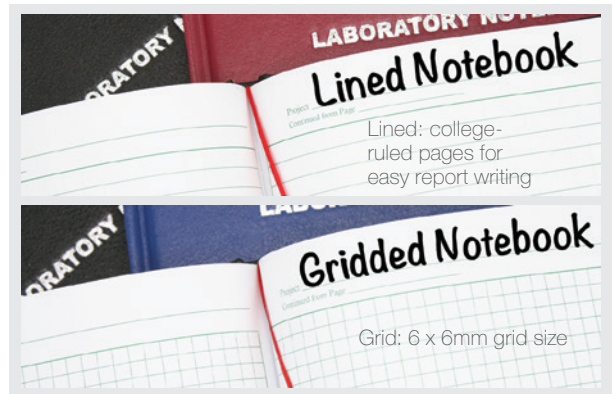

Fisherbrand™ Laboratory Notebooks

Keep information safe and well organised. Record research data, experimental procedures, design and development ideas, and more.

- Smyth Sewn binding for security; pages lie flat when notebook is open
- Acid-free paper pages are numbered, with spaces for title, date, project number, book number, date, and witness signature
- Each notebook includes instructions, table of contents, periodic table, and conversion tables

| Cat. No. | Colour | Pages | Page Format. | L x W cm | Qty |
|----------|--------|-------|--------------|-------------|-----|
| 15801658 | Black | 100 | Grid | 21.8 x 28.1 | 6 |
| 15841658 | Blue | | | | |
| 15881658 | Red | | | | |
| 16216870 | Green | 200 | Lined | | |
| 15811658 | Black | | | | |
| 15851658 | Blue | | | | |
| 15891658 | Red | 100 | Lined | | |
| 16226870 | Green | | | | |
| 15821658 | Black | | | | |
| 15861658 | Blue | 200 | Lined | | |
| 15801668 | Red | | | | |
| 16236870 | Green | | | | |
| 15831658 | Black | 100 | Lined | | |
| 15871658 | Blue | | | | |
| 15811668 | Red | | | | |
| 16246870 | Green | 200 | Lined | | |
| | Black | | | | |
| | Blue | | | | |

NEW

Fisherbrand™ Lab Drying Rack

Bringing colour to a mundane task

Most versatile rack in the market


- Channels allow for more useful configurations
- Large and small pegs for all glassware varieties
- Unique drying basket to collect smaller items

Designed to protect against breakage

- Angled peg allows for better hold - especially large, heavy pieces
- Patent pending SureGrip hook securely holds glassware in place

Included:

- 72 Pegs: 36 small/green, 36 large/blue
- 2 drying baskets
- Drainage tubing
- Mounting hardware

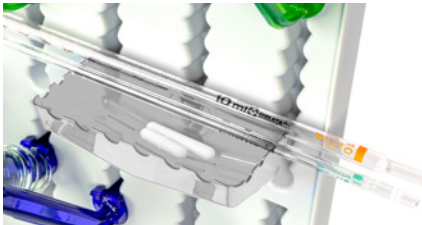

Drying Basket

Keeps small items and pipets safe



SureGrip Hook

Patent-pending design securely grips glassware into place




Use on counter or wall

Stands alone on integrated hinged stand. Mounts easily on wall with included guide and hardware

Can hold up to 34.9 kg. (77 lbs)

Rack made of durable polystyrene

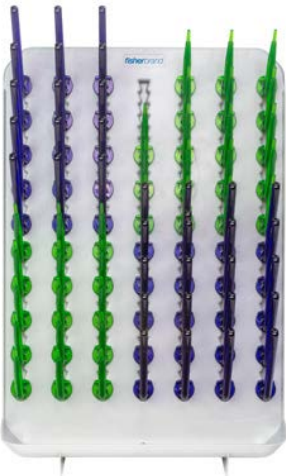

Lab Drying Rack

| Cat. No. | L x W x H cm | Qty |
|----------|------------------|-----|
| 17120041 | 44.1 x 66 x 13.9 | 1 |

Replacement Parts

| Cat. No. | Description | Qty |
|----------|------------------------------------------|-----|
| 17190031 | Peg 10 Piece Kit, 5 Small, 5 Large | 10 |
| 17130041 | Drying Basket Kit | 2 |
| 17100041 | Replacement Draining Set, Tube and Elbow | 1 |

Fisherbrand™ Wire Drying Rack

Saves glassware from chips and cracks

Drain and dry test tubes, beakers, flasks and other labware

- Rack is self-standing; or use integral hooks to hang on the wall
- Holds up to 90 items on 18 x 127 mm pins and 72 x 76 mm pins
- Made from HDPE-coated spot-welded steel wire for durability; polystyrene back
- Easily accessible and removable polystyrene drip pan



Removable Drip Pan
for water collection and disposal

| Cat. No. | L x W x H cm | Qty |
|----------|--------------------|-----|
| 16330746 | 46.2 x 18.2 x 52.5 | 1 |

Fisherbrand™ Safety Bottle Carrier

Secure and convenient bottle transport

- Ergonomic handle prevents hand fatigue and related accidents when carrying heavier loads
- Handle swings down to either side for easier loading and unloading
- Centering ring stabilizes bottle in the center of carrier for safer, steadier transport
- Carriers come in four colors with designated labeling areas for easy identification
- Made of chemical resistant thermoplastic elastomer - PP blend

NEW



Available in 4 Colors
for easy identification and organization

Safety Bottle Carrier, 4.5 L


| Cat. No. | Color | Dimensions (Dia. x H) in | Dimensions (Dia. x H) cm | Qty |
|----------|-------|-----------------------------------------------|--------------------------------------------------|---------------|
| 17840782 | Red | 7.2 x 10.5 (16 in height, handle extended) | 18.3 x 26.7 (40.6 cm height, handle extended) | 1 ea, 4 ea/cs |
| 17830782 | Black | | | |
| 17850782 | Blue | | | |
| 17810782 | Green | | | |

Safety Bottle Carrier, 1.5 L

| Cat. No. | Color | Dimensions (Dia. x H) in | Dimensions (Dia. x H) cm | Qty |
|----------|-------|----------------------------------------------|--------------------------------------------------|---------------|
| 17820782 | Red | 4.5 x 6.5 (10 in height, handle extended) | 11.4 x 16.5 (25.4 cm height, handle extended) | 1 ea, 4 ea/cs |
| 17800782 | Black | | | |
| 17890772 | Blue | | | |
| 17880772 | Green | | | |



Movable Handle
swings down to allow clear access to vessel

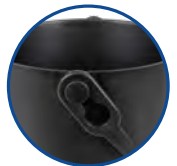



Fisherbrand™ RF1000

Lightweight and easy to use

The RF1000 manual pipette controller is uniquely designed for one-handed operation

Simply squeeze the main bulb and push the toggle forward to aspirate or backward to dispense. Push the blowout bulb to expel any liquid remaining in the pipette.

- Tapered silicone pipette adapter works with 1 to 100 mL pipets and is replaceable
- Includes a replaceable 0.45 µm hydrophobic filter to protect your samples and internal components
- One-year warranty
- **Autoclavable** nose cone assembly only

RF1000

| Cat. No. | Colour | Volume Range mL | L x W x H cm | Qty |
|----------|--------|-----------------|---------------|-----|
| 11821920 | Blue ■ | 1-100 | 18 x 4.5 x 12 | 1 |

Replacement Parts

| Cat. No. | Description | Material | Colour | Qty |
|----------|-----------------------------|----------|--------|-----|
| 11304035 | Replacement Filters | | | 5 |
| 11831920 | Replacement Pipette Adapter | Silicone | | 1 |
| 11891920 | Replacement Nose Cone | PP | Blue ■ | 1 |



Stay Organised

with two additional nose cones (green, white) for colour identification for both the Pipette Controller & RF1000

Fisherbrand™ Fillers And Bulbs

Rubber Bulbs (Natural Rubber)


| Cat. No. | Colour | Volume mL | Length mm | Outer Diameter mm | Hole Diameter mm | Qty |
|----------|--------|-----------|-----------|-------------------|------------------|-----|
| 11977984 | Red ■ | 15 | 56 | 34 | 4 | 1 |
| 11987984 | | 30 | 66 | 46 | 6 | 1 |
| 11997984 | | 60 | 80 | 54 | 6 | 1 |
| 11314055 | | 150 | 103 | 73 | 6 | 1 |

Contoured shape and heavy-wall neck increase grip for easier use with pipe

Safety Bulb Pipette Filler (Natural Rubber)

| Cat. No. | Colour | Volume mL | Length mm | Outer Diameter mm | Hole Diameter mm | Qty |
|----------|--------|-----------|-----------|-------------------|------------------|-----|
| 11779165 | Red ■ | 15 | 106 | 41 | 4.5-6.5 | 3 |

- Features a white ABS attachment containing a silicone conical adapter that seals around pipettes without twisting
- Pipette does not contact hard plastic, which reduces the risk of cracking and chipping
- Elastic band for hanging on reagent bottle or wrist
- Filler fits standard size pipettes

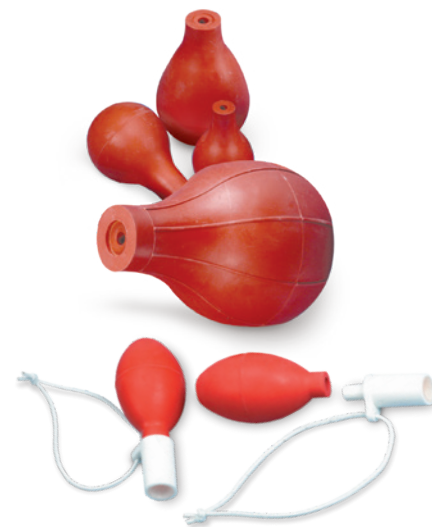

Fisherbrand™ Reservoir Chillers

Modern Design & Enhanced Features

Ice-free, non-toxic way to extend working period and preserve integrity for temperature sensitive reagents

- Keeps reagents at temperatures less than 5°C (41°F) for approximately 3 hours
- Engineered custom fit design for 50/55 or 100 mL reservoirs for maximum contact and for longer uniform cooling
- Eliminates the risk of splash contamination from ice baths

Stackable
Saves storage and work space



NEW

Symmetrical cooler design

Allows reagent reservoirs to be mounted in either direction

50/55 mL



100 mL



| Cat. No. | Colour | Volume mL | L x W x H cm | Qty |
|----------|-----------|-----------|------------------|-----|
| 16330412 | Purple | 50/55 | 14.4 x 7.4 x 4.5 | 2 |
| 16340412 | Blue | | | |
| 16350412 | Assorted* | 100 | 15 x 9 x 4.3 | 2 |
| 16360412 | Purple | | | |
| 16370412 | Blue | | | |
| 16380412 | Assorted* | | | |

*Assorted colours (Purple and Blue)

Fisherbrand™ Multi-Channel Reagent Reservoirs

Designed to streamline pipetting procedures

- Superior durability with enhanced walls and crush-resistant "turrets" at the corners
- Easy handling with tabs to prevent stacked reservoirs from sticking together
- Sterile versions are gamma irradiated in an ISO-certified facility and individually sealed in HDPE packaging; traceability certificate available upon request

| Cat. No. | Volume mL | Material | Sterility | L x W x H cm | Qty |
|----------|-----------|-------------------|-------------|---------------------|-----|
| 16371076 | 55 mL | Clear PVC | Non-Sterile | 10.2 x 2.8 x 3.0 cm | 100 |
| 16361076 | 55 mL | White Polystyrene | Non-Sterile | | |
| 16351076 | 55 mL | White Polystyrene | Sterile | 10.9 x 4.8 x 3.2 cm | 100 |
| 16344862 | 100 mL | White Polystyrene | Non-Sterile | | |
| 16334862 | 100 mL | White Polystyrene | Sterile | | |

55 mL




100 mL



Easy pour spouts

Tapered "V" Design

Allows for full sample recovery

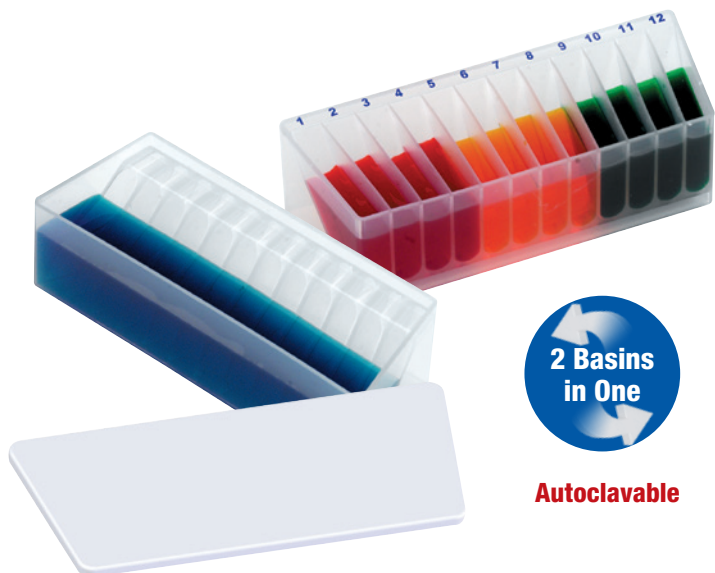

Autoclavable

Fisherbrand™ Dual Solution Reservoir

Innovative shape guides and protects pipette tips

Multi-functional reagent reservoir has a single 50 mL basin on one side and 12 slotted 5 mL basins on the other side

- Angled wells guide pipette tips toward a centre position to help avoid cross contamination and waste
- Individual 5 mL basin wells numbered for sample identification
- Covers (included) help prevent contamination

| Cat. No. | Volume mL | Sterility | L x W x H cm | Qty |
|----------|-----------|-------------|----------------|-----|
| 11978084 | 5/50 | Non-Sterile | 12.5 x 5 x 3.7 | 25 |

Fisherbrand™ Universal Pipette Stand

Pipettes never looked so good on your bench top. Stylish pipette stand holds and stores up to 8 single-channel, 4 multichannel, or a combination of pipettes.

Versatile

- Patent pending design securely holds virtually every brand of pipette in the marketplace
- Modular upper holder inserts accommodate different pipette widths. Inserts come in two separate colours allowing for colour coding configurations

Durable Utility

- Stable weighted base allows for easy one-handed pipette removal and replacement
- Rotating head makes pipette selection fast and easy
- Quick grip handle lets you move the stand while keeping the pipettes in place





Made of Durable ABS

Easy to wipe clean



Vertical Storage

Helps prevent contamination

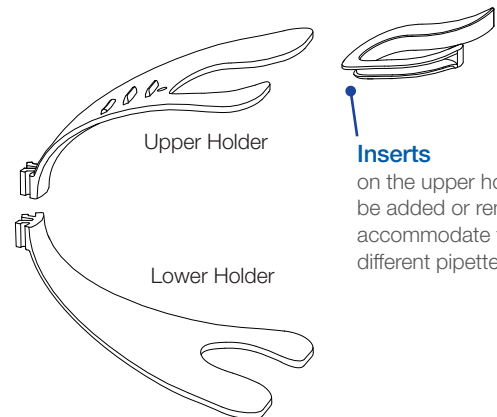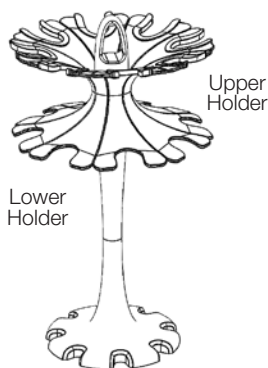

Inserts

on the upper holders can be added or removed to accommodate the fit of different pipette widths

US Patent Pending

| Cat. No. | Upper Holder Colour | Quantity | Lower Holder Colour | Quantity | Base Width cm | Top W x H cm | Qty |
|----------|-------------------------------------------------------------------------------------|---------------|------------------------------------------|----------|---------------|--------------|-----|
| 15650785 | Blue/Green ■ ■ | 8 each/colour | Blue ■ | 8 each | 15.8 | 22.9 x 40.8 | 1 |

Fisherbrand™ Universal Pipette Wall-Mount

Saves valuable bench top space

The pipette holder comes with two wall mounting brackets that each hold two pipettes.

Secure

Upper and lower segment inserts easily slide in place to securely hold virtually every brand of pipette


Unique design includes traction grips that centre pipettes for one-handed removal and replacement

Adaptable

- Inserts come in two separate colours allowing for colour coding and quick identification or add a pop of colour to your workspace



Vertical Storage
Helps prevent contamination

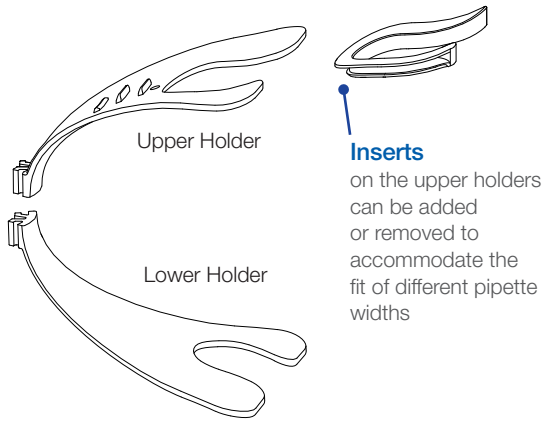

US Patent Pending

| Cat. No. | Lower Holder | Upper Holder | L x H cm | Qty |
|----------|--------------|--------------|------------|-----|
| 15847944 | 2 Blue | 2 Blue | 17.5 x 4.5 | 2 |
| | 2 Green | 2 Green | | |
| | 2 Clear | 1 Purple | | |
| | | 1 Yellow | | |



Fisherbrand™ Pipette Holders

Secure, easy clamping

- Triple format fits most major pipette brands
- Made of chemical-resistant polypropylene

Autoclavable

| Cat. No. | Colour | L x W x H cm | Qty |
|----------|--------|-----------------|-----|
| 11821582 | Blue | 15.7 x 12 x 9.5 | 1 |



Secure Clamping
for tabletops, benches and shelves



Fisherbrand™ Pipette Holder

Giving labs a choice

Unique pipette holder gives you the option of using on the bench top, or hanging under a shelf to help maximise your working space

- Holds 12 single-channel, or a combination of 3 single-channel and 3 multichannel pipettes
- Weighted base and rotating rack make for easy pipette selection
- Rubber grips prevent pipette damage

| Cat. No. | Colour | Bench Top Style W x H cm | Under Shelf Style W x H cm | Qty |
|----------|---------|--------------------------|----------------------------|--------|
| 15291588 | Blue | | | |
| 15380360 | Lime | 33.2 x 19.3 | 18.7 x 19.3 | 1 each |
| 15390360 | Magenta | | | |



Pipettes can be held in the rubber grips and/or the plunger hooks



Mount base undershelf or counter with included hardware



US Design Patent No. D754872
US Design Patent No. D763464



3-Place Stand Acrylic Clear



4-Place Stand Acrylic Clear



6-Place Stand Acrylic Blue



6-Place Stand ABS Blue

Fisherbrand™ Pipette Stands

Convenient and safe storage

Pipette stands are designed to hold most major pipette brands

- Acrylic stands are transparent/clear
- ABS stands are opaque

| Cat. No. | Colour | Material | Pipette Capacity | Slot Upper mm | Slot Lower mm | Qty |
|----------|------------------------------------------|----------|------------------|---------------|---------------|-----|
| 10760774 | Clear <input type="checkbox"/> | Acrylic | 3 | | | |
| 10122952 | Clear <input type="checkbox"/> | Acrylic | 4 | | | |
| 11394045 | Clear <input type="checkbox"/> | Acrylic | 6 | 25 | 20 | 1 |
| 11384045 | Blue <input checked="" type="checkbox"/> | Acrylic | 6 | | | |
| 10563893 | Blue <input checked="" type="checkbox"/> | ABS | 6 | | | |

Fisherbrand™ Cargo Pipet Rack

Meets the needs of any lab

Customisable for any situation

Racks securely lock together horizontally or vertically or can be used individually

Place on a benchtop or shelf or use the built-in magnets to attach to a steel surface


- Angled shelves store 0.1 to 50 mL pipettes
- Protective lid encloses the unit to prevent dust build-up and allows stacking

Made of Durable
ABS
Easy to wipe clean



| Cat. No. | Colour | L x W x H cm | Qty |
|----------|--------|------------------|-----|
| 11958004 | Blue | 30.1 x 10 x 12.8 | 3 |

Magnetic
Attaches to any
steel surface with
strong built-in
magnets

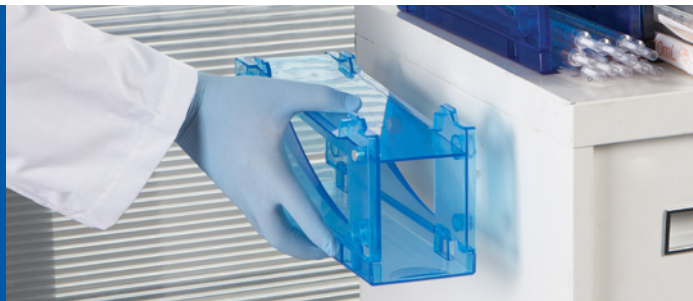


Fisherbrand™ Pipet Racks

Protects and organises glass and plastic pipettes

Four angled-shelved compartments keep your 0.1 to 50 mL pipets or cans within easy reach while preventing damage and clutter

- Heavy-duty acrylic version provides visibility from any angle
- ABS versions have cut-out viewing windows and are supplied with two removable holders for bulbs, fillers, and pumps
- HS20612M magnetised ABS version mount securely to steel surfaces with magnets on both sides

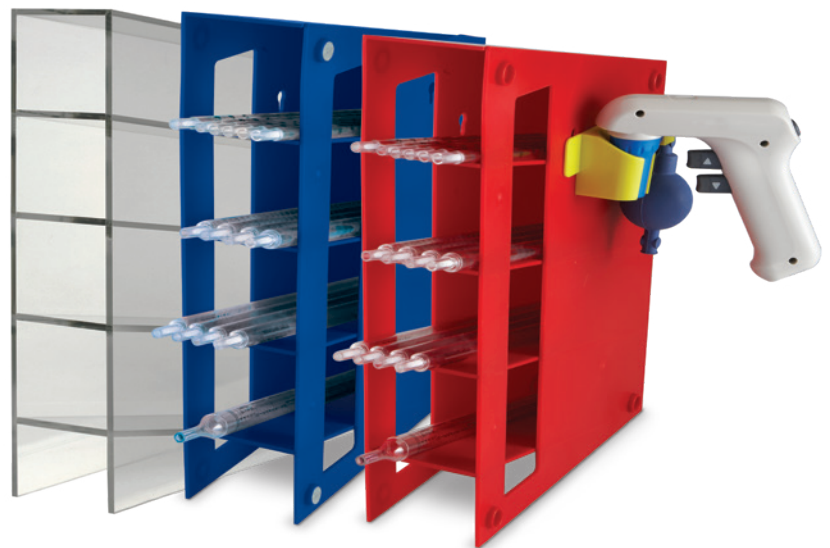

| Cat. No. | Colour | Material | L x W x H cm | Qty |
|----------|--------|------------------------|-------------------|-----|
| 11703926 | Clear | Acrylic | | |
| 11793849 | Blue | ABS | 29.1 x 8.8 x 40.1 | 1 |
| 11791455 | Red | ABS | | |
| 11757325 | Blue | ABS (Mounting Magnets) | | |



Pipette Attachments
ABS version includes 2
removable attachments

NEW

Fisherbrand™ Clinical 4-Way Tube Rack

Multiple combinations for any situation

These racks can be connected from all six sides to create customised configurations and hold round-, flat-, or conical-bottomed tubes.

Versatile

- Use as standalone unit or link multiple units to fit your specific needs
- Holds any tube with a diameter between 13 and 17 mm
- Assorted pack contains five differently coloured racks

Holds Multiple Tube Sizes without multiple racks




Robust Functional Design

- Uniform well spacing helps eliminate pipetting errors for more accurate results
- Moulded alphanumeric identification for easy tube location and retrieval
- Made from chemically resistant and easy-to-clean polypropylene
- Racks can be used in water baths or placed in freezers

Easy to Configure

Racks slide together to create many combinations

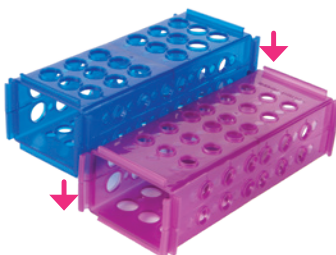

Autoclavable

Tube Bottoms Nest

in well pockets keeping them from going thru the rack and off the bench top


Holds 15 mL Conical Tubes



Patent Pending

| Cat. No. | Colour | Cat. No. | Colour |
|----------|--------|----------|-----------|
| 16350402 | Blue | 16380402 | Natural |
| 16360402 | Purple | 16340402 | Assorted* |
| 16370402 | Green | | |

*Assorted colours (Blue, Purple, Green, Red, Natural)

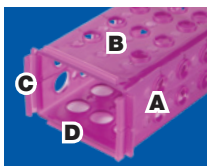

Individual Rack Specifications

| Side | Wells | Well Dia. mm | L x W x H cm | Qty |
|------|-------|--------------|------------------|-----|
| A | 12 | 13 | 23.2 x 7.0 x 9.5 | 5 |
| B | 18 | 13 | | |
| C | 12 | 17 | | |
| D | 18 | 17 | | |



Fisherbrand™ 4-Way Tube Rack

Link together to meet your application needs

Versatile

- Use as standalone unit or link multiple units to fit your specific needs
- Holds tubes with diameters from 8.2 to 30 mm


Robust Functional Design

- Moulded alphanumeric identification for easy tube location and retrieval
- Made from chemical resistant and easy to clean polypropylene
- Racks can be used in water baths or placed in freezers

Easy to Configure

Racks slide together end to end or side to side to create many combinations

Autoclavable

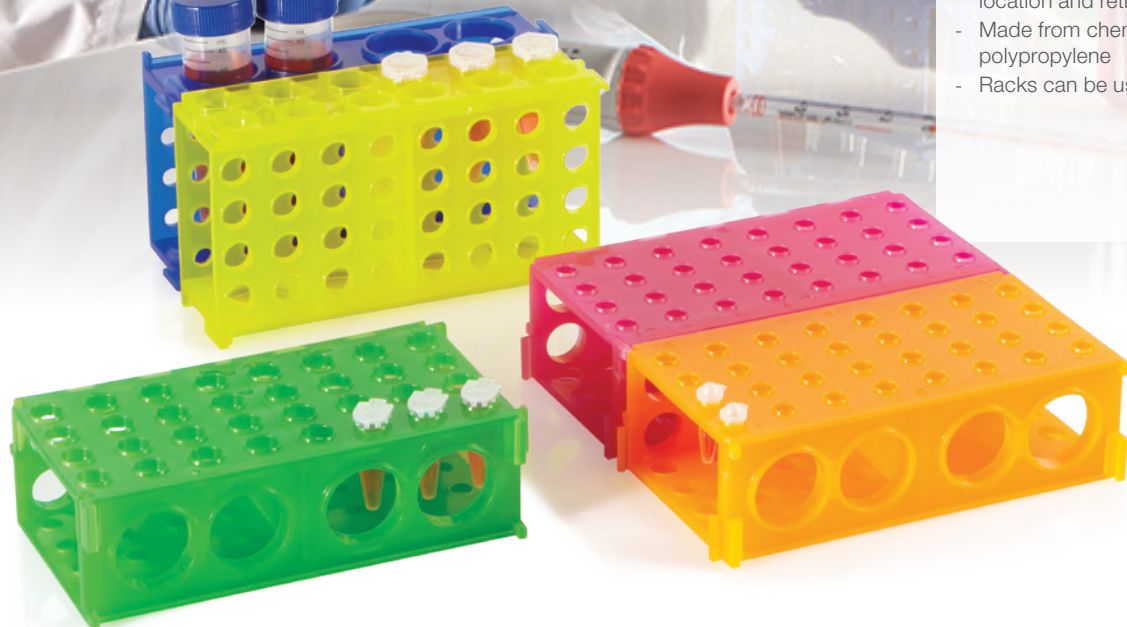

| Cat. No. | Colour | Cat. No. | Colour | Tubes mL | Wells | Well Dia. mm | L x W x H cm | Qty |
|----------|---------|-------------------------------------------------------|-----------|----------|---------|--------------|--------------------|-----|
| 11750334 | Natural | 11770334 | Yellow | 50 | 4 | 30 | 10.8 x 10.8 x 11.3 | 5 |
| 11770045 | Blue | 11790045 | Orange | 15 | 12 | 17.5 | | |
| 11760334 | Green | 11700055 | Assorted* | 0.5 | 32 | 8.2 | | |
| 11780045 | Pink | *Assorted Colours (Blue, Green, Pink, Yellow, Orange) | | | 1.5/2.0 | 32 | | |

Fisherbrand™ Cube Tube Rack

Interlocking multi-tube racks

Create custom configurations can be made by linking cubes together

- Assorted colours make colour coding easy
- Made from durable polypropylene that withstands ice baths and freezer conditions and can be easily cleaned



| Cat. No. | Colour | Tubes mL | Wells | Well Dia. mm | L x W x H cm | Qty |
|----------|--------|----------|-------|--------------|--------------------|-----|
| 11784726 | | 50 | 4 | 30 | 10.8 x 10.8 x 11.3 | 5 |
| | | 15 | 10 | 17 | | |
| | | 3 to 7 | 12 | 12.5 | | |
| | | 1.5/2.0 | 16 | 11.5 | | |

*Assorted colours (Blue, Green, Pink, Yellow, Orange)

Autoclavable

Fisherbrand™ Multi Tube Rack

Exclusive design, optimal functionality

This fully flexible tube rack is designed for tubes of any height with diameters between 10 and 30 mm.


Multiple Configurations for Workflow Efficiency

- Interchangeable tube rack inserts hold multiple tube sizes firmly in place when alternating between sample vessels
- Inserts can be pivoted in either an upper or lower position to securely accommodate the size and height of tube being used

Multi-Functional Design

- Open design allows for better filling visibility of tube graduations when filling
- Rack wells have a unique teardrop design that enables better retention of tubes
- Rotate tube rack insert any angle for more comfortable pipetting

Unique Well Design
Tear drop shape holds tube snugly

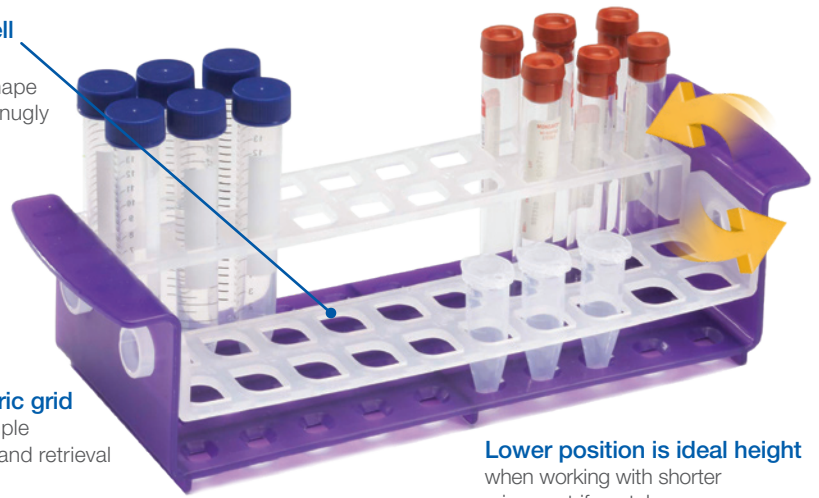

Alphanumeric grid for quick sample identification and retrieval

Lower position is ideal height when working with shorter microcentrifuge tubes

Ergonomic handles for ease of transport



Angled position for more comfortable pipetting




| Included: | Two 10 to 13 mm Tube Racks Fits | Two 13 to 17 mm Tube Racks Fits | Two 25 to 30 mm Tube Racks Fits |
|-------------------------------------|---------------------------------|---------------------------------|---------------------------------|
| 1 Base + 6 Tube Rack Inserts | 5 mL, 7.5 mL, 10 mL Tubes | 5 mL, 15 mL Tubes | 25 mL, 50 mL Tubes |

Multi Tube Rack

US Patent No. 9776188

| Cat. No. | Colour | L x W x H cm | # Inserts Included | Well Dia./Insert | Wells/Insert | Qty |
|----------|----------------|-------------------|--------------------|------------------|--------------|-----|
| 17181902 | Blue/Natural | 30.9 x 10.7 x 7.7 | 6 (2 ea/ Diameter) | 13 mm | 20 | 1 |
| | | | | 17 mm | 20 | |
| | | | | 30 mm | 6 | |
| 17191902 | Purple/Natural | 30.9 x 10.7 x 7.7 | 6 (2 ea/ Diameter) | 13 mm | 20 | 1 |
| | | | | 17 mm | 20 | |
| | | | | 30 mm | 6 | |

Alphanumeric well identification

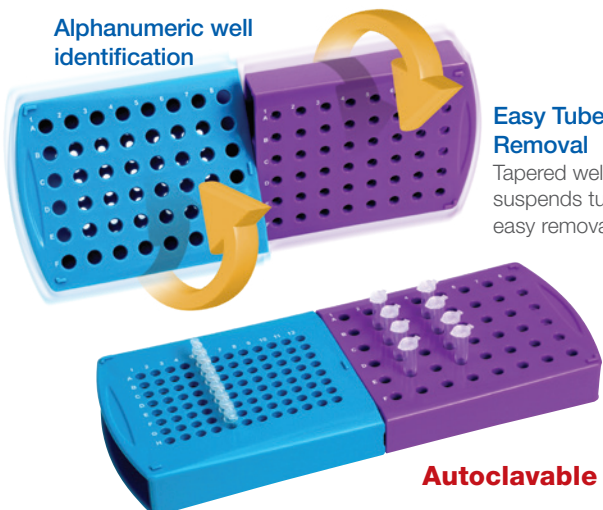

Easy Tube Removal
Tapered well design suspends tube for easy removal

Autoclavable

Fisherbrand™ Duo Tube Rack

Rotate to customise workflow processing

Transfer solutions between 0.2 mL PCR® tubes and 0.5 and 1.5 mL microcentrifuge tubes

- Purple section holds 48 x 0.5 mL and 48 x 1.5/2.0 mL microfuge tubes on the opposite side
- Blue section holds 96 x 0.2 mL tubes and 48 x 1.5/2.0 mL tubes on the opposite side

| Cat. No. | Colour | Tubes mL | Wells | Well Dia. mm | L x W x H cm | Qty |
|----------|-------------|----------|-------|--------------|-------------------|-----|
| 11304095 | Blue/Purple | 0.2 | 96 | 5 | 29.6 x 11.8 x 3.4 | 1 |
| | | 0.5 | 48 | 8 | | |
| | | 1.5/2.0 | 96 | 11 | | |


Fisherbrand™ Tube Racks

Modular, fun rotating rack

Each module is interlocked and adjustable with just a twist allowing you to work with multiple sized tubes in the same rack

- Suitable for variety of culture, centrifuge, micro-centrifuge, round, flat or conical bottomed tubes
- Two versions to accommodate small and large tube sizes

Alphanumeric well identification

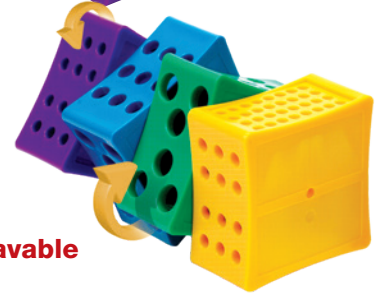

Autoclavable

| Cat. No. | Colour | Tubes mL | Wells | Well Dia. mm | L x W x H cm | Qty |
|----------|--------|----------|-------|--------------|----------------|-----|
| 11384085 | Yellow | 50 | 2 | 29 | 22.4 x 11 x 11 | 1 |
| | Green | 10 to 15 | 8 | 20 | | |
| | Blue | 15 to 20 | 6 | 16 | | |
| | Purple | 10 to 13 | 10 | 13.5 | | |
| 11394085 | Yellow | 15 | 6 | 16.4 | 21 x 10 x 10 | 1 |
| | Green | 1.5/2.0 | 9 | 11.2 | | |
| | Blue | 0.5 | 12 | 8 | | |
| | Purple | 0.2 | 32 | 6.5 | | |

Fisherbrand™ Adaptable Multi-tube Rack

Flexible & linkable to your needs

Made from POM acetyl plastic with flexible fingers that adapt to accommodate 12 to 30 mm tube diameters and hold them firmly in place.

- Racks can be linked together to expand capacity
- Open sides allow for better label and fill level viewing
- Available in single or dual colours for easy identification

**Fits
25 mL
Conical
Tubes**



Autoclavable

| Cat. No. | Rack #1 Colour | Rack #2 Colour | Tubes mL | Wells | Well Dia. mm | L x W X H cm | Qty |
|----------|----------------|----------------|----------|-------|--------------|------------------|-----|
| 15300370 | Blue/Green | Green/Blue | 7 to 50 | 4 | 12 to 30 | 18.1 x 5.6 x 7.6 | 2 |
| 15320370 | Blue/Yellow | Yellow/Blue | | | | | |
| 15330370 | Pink/White | White/Pink | | | | | |
| 15340370 | Blue/Blue | Blue/Blue | | | | | |
| 15350370 | White/White | White/White | | | | | |
| 15360370 | Pink/Pink | Pink/Pink | | | | | |
| 15380370 | Yellow/Yellow | Yellow/Yellow | | | | | |
| 15390370 | Green/Green | Green/Green | | | | | |

Fisherbrand™ One-Well Connecting Tube Racks

Holds 50 mL tubes

Form a line or use as a single, separate rack

- Open design enables easy sample identification
- Moulded circular base ring keeps conical or round-bottom tubes in place

| Cat. No. | Colour | Tubes mL | Wells | Well Dia. mm | L x W x H cm | Qty |
|----------|------------------------------------------------------|----------|-------|--------------|---------------|-----|
| 11727285 | Assorted colours (Blue, Green, Yellow, Orange, Pink) | 50 | 1 | 29 | 8.3 x 7 x 6.7 | 5 |

*Assorted colours (Blue, Green, Yellow, Orange, Pink)




Fisherbrand™ Two Tier Microtube Rack 20-Well

Autoclavable

Interlock multiple racks to increase tube storage

- Tiered rack allows for improved staging and visibility for easier, more accurate sample filling and retrieval
- Made of easy-to-clean polypropylene
- Easily remove tubes using finger slot on front of the rack
- Embossed index numbers for quick sample ID and retrieval

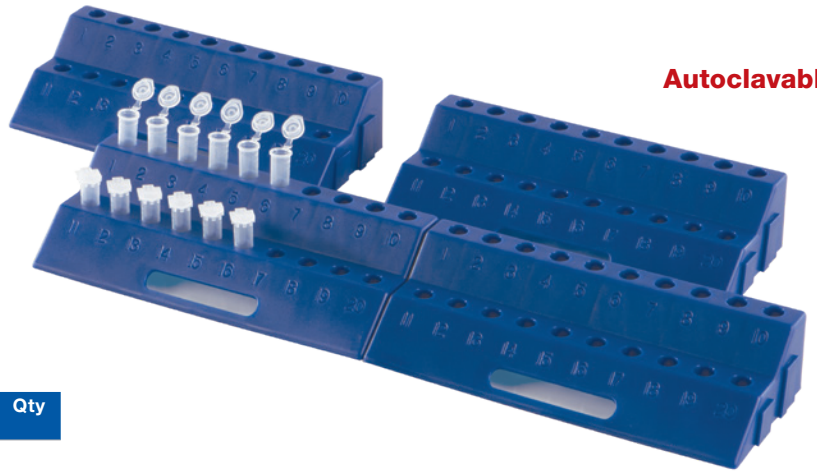


| Cat. No. | Colour | Tubes mL | Wells | Well Dia. mm | L x W x H cm | Qty |
|----------|--------|----------|-------|--------------|--------------|-----|
| 11784016 | Blue | 1.0 | 20 | 10 | 21.1 x 9 x 5 | 6 |
| | | 1.5 | | | | |
| | | 2.0 | | | | |

Fisherbrand™ Universal Tube Rack

Ideal for holding a wide range of tubes

- Holds up to 12 similar size tubes
- Alphanumeric identification for quick location and retrieval
- Side handles for easy transport
- Made from polypropylene

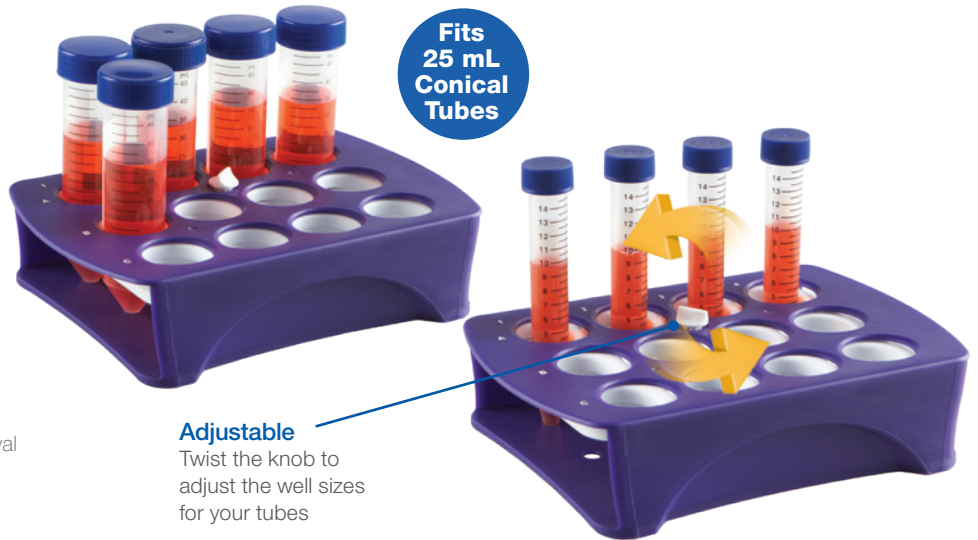

| Cat. No. | Colour | Tubes mL | Wells | Well Dia. mm | L x W x H cm | Qty |
|----------|--------------|----------|-------|--------------|-----------------|-----|
| 11987994 | White/Purple | 15 to 50 | 12 | 16 to 30 | 19.4 x 14 x 5.1 | 1 |

Fisherbrand™ Transformable Cube Tube Rack

Autoclavable

Lidded tube rack & storage

Quick change over with 2 snap-in inserts that hold 15 and 50 mL tubes

- Keyed lid and numeric wells for tube identification
- Made from polypropylene to withstand storage temperatures down to -80°C

Compact Storage

- Clear Lid included for viewable, stackable storage
- Space available under insert for ice or extra storage

| Cat. No. | Colour | Tubes mL | Wells | Well Dia mm | L x W x H cm | Qty |
|----------|--------|----------|-------|-------------|--------------------|-----|
| 11700634 | Purple | 15 | 16 | 18 | 12.5 x 12.5 x 12.9 | 5 |
| | | 50 | 9 | 31 | | |





Fisherbrand™ Snap-Together Conical Tube Rack

Space Saver - shipped flat and easy-to-assemble

Conveniently holds 15 and 50 mL tubes simultaneously

- ABS racks snap together securely and are stable in waterbaths
- Alphanumeric identification for quick location and retrieval

| Cat. No. | Colour | Tubes mL | Wells | Well Dia. mm | L x W x H cm | Qty |
|----------|--------|----------|-------|--------------|-----------------|-----|
| 10753257 | Orange | 15 | 30 | 17 | 17.1 x 17 x 6.1 | 1 |
| | | 50 | 20 | 28 | | |

Fisherbrand™ Pop-Up Tube Racks

Collapse for compact storage and “pop-up” for quick use

Designed to save space

- Easy to open, holds 15 and 50 mL tubes
- Easy to close, stackable for compact storage
- Alphanumeric well identification for quick location and retrieval of samples



| Cat. No. | Colour | Tubes mL | Wells | Well Dia. mm |
|----------|--------|----------|--------|--------------|
| 11710634 | Blue | 15, 50 | 21, 12 | 17, 30 |
| 12269309 | Pink | | | |
| 12693466 | Purple | 15 | 45 | 17 |
| 12603476 | Green | 50 | 18 | 30 |

| L x W x H (Expanded) cm | L x W x H (Collapsed) cm | Qty |
|-------------------------|--------------------------|-----|
| 25.5 x 13.7 x 7.2 | 25.5 x 13.7 x 2.1 | 2 |



Fisherbrand™ Tube Racks

Tough and highly resistant to acids, bases, solvents and heat

Durable racks are moulded in a single, continuous piece so no assembly required and no detachable pieces

Intuitive Design

- Unique tear drop shaped openings offer better tube retention
- Moulded in coordinates for tube orientation
- Flat labelling area for easier identification
- Stable design for use with water baths

Available in two materials:


Fiberglass-Reinforced Polyoxymethylene (POM)

- Available in half and full rack sizes for 5 different tube diameters
- Offers the highest chemical and heat resistance (-40 to +121°C)

Autoclavable

Polypropylene

- Available in full rack sizes for 2 different tube diameters
- Offers good chemical and heat resistance (0 to 110°C)

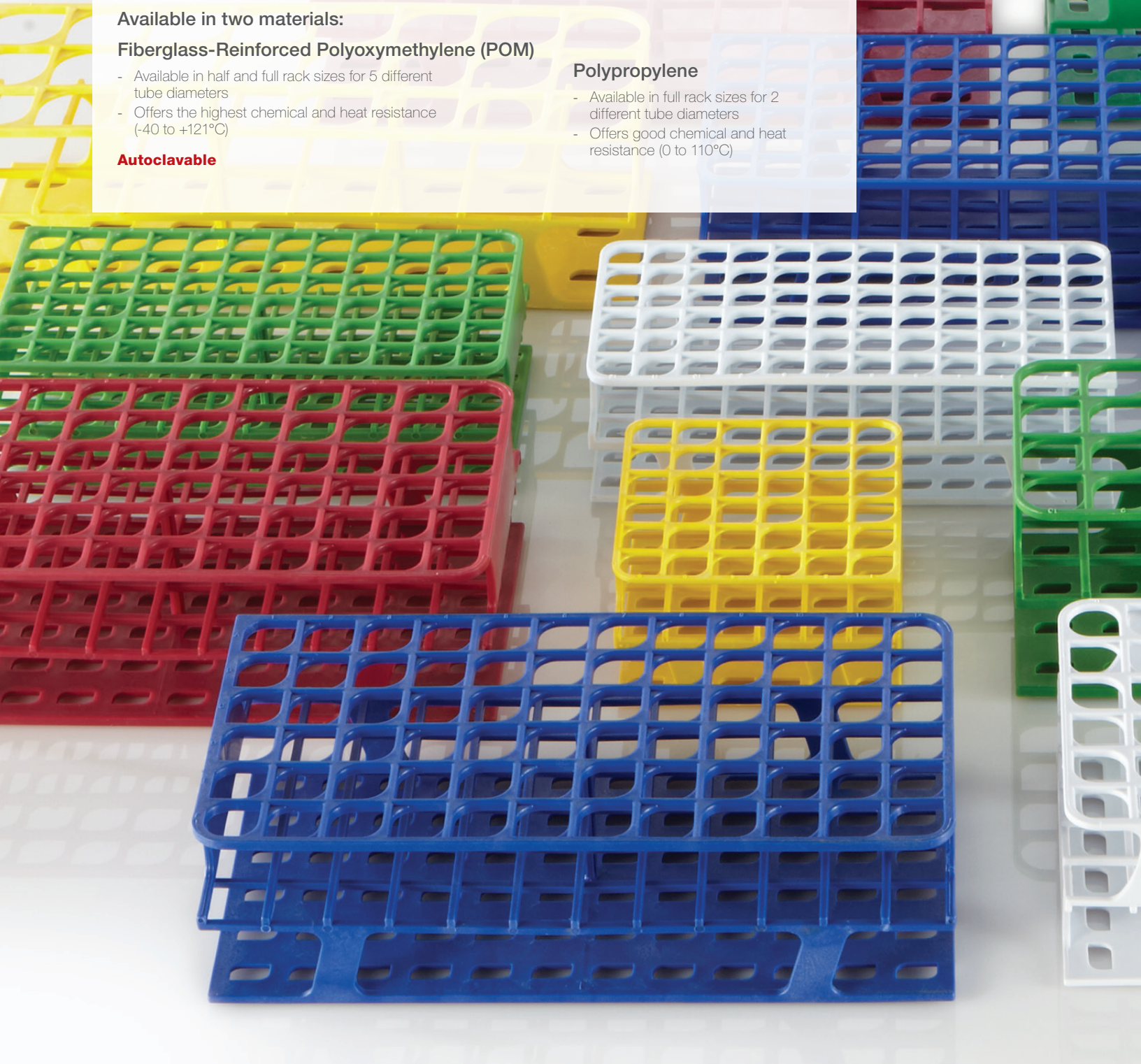



Full-Size POM

| Cat. No. | Colour | Tubes mL | Wells | Well Dia. mm | L x W x H cm | Qty |
|----------|----------------------------------|----------|-------|--------------|-------------------|-----|
| 15340380 | White <input type="checkbox"/> | 5 to 10 | 72 | 13 | 20.2 x 10.4 x 5.9 | 1 |
| 15350380 | Blue <input type="checkbox"/> | | | | | |
| 15360380 | Green <input type="checkbox"/> | | | | | |
| 15370380 | Magenta <input type="checkbox"/> | | | | | |
| 15380380 | Yellow <input type="checkbox"/> | | | | | |
| 15390380 | White <input type="checkbox"/> | 5 to 10 | 72 | 16 | 25 x 12.7 x 7 | 1 |
| 15300390 | Blue <input type="checkbox"/> | | | | | |
| 15310390 | Green <input type="checkbox"/> | | | | | |
| 15320390 | Magenta <input type="checkbox"/> | | | | | |
| 15330390 | Yellow <input type="checkbox"/> | | | | | |
| 15340390 | White <input type="checkbox"/> | 10 to 18 | 40 | 20 | 25.2 x 10 x 8.3 | 1 |
| 12327629 | Blue <input type="checkbox"/> | | | | | |
| 15350390 | Green <input type="checkbox"/> | | | | | |
| 15360390 | Magenta <input type="checkbox"/> | | | | | |
| 15370390 | White <input type="checkbox"/> | | | | | |
| 15380390 | Blue <input type="checkbox"/> | 10 to 18 | 40 | 25 | 30 x 12 x 9.2 | 1 |
| 15390390 | Green <input type="checkbox"/> | | | | | |
| 15300400 | Magenta <input type="checkbox"/> | | | | | |
| 15310400 | White <input type="checkbox"/> | | | | | |
| 10257963 | Blue <input type="checkbox"/> | | | | | |
| 15330400 | Green <input type="checkbox"/> | 15 to 50 | 24 | 30 | 28.2 x 11 x 8.5 | 1 |
| 15350400 | Magenta <input type="checkbox"/> | | | | | |
| 15360400 | Yellow <input type="checkbox"/> | | | | | |

Half-Size POM

| Cat. No. | Colour | Tubes mL | Wells | Well Dia. mm | L x W x H cm | Qty |
|----------|----------------------------------|----------|-------|--------------|-------------------|-----|
| 15390400 | White <input type="checkbox"/> | 5 to 10 | 36 | 13 | 10.4 x 10.4 x 5.9 | 1 |
| 15300410 | Blue <input type="checkbox"/> | | | | | |
| 15310410 | Green <input type="checkbox"/> | | | | | |
| 15330410 | Magenta <input type="checkbox"/> | | | | | |
| 15519821 | Yellow <input type="checkbox"/> | | | | | |
| 15340410 | White <input type="checkbox"/> | 5 to 10 | 36 | 16 | 12.7 x 12.7 x 7 | 1 |
| 15360410 | Blue <input type="checkbox"/> | | | | | |
| 15380410 | Green <input type="checkbox"/> | | | | | |
| 15390410 | Magenta <input type="checkbox"/> | | | | | |
| 15529821 | Yellow <input type="checkbox"/> | | | | | |
| 15310420 | White <input type="checkbox"/> | 10 to 18 | 20 | 20 | 10 x 12.7 x 8.3 | 1 |
| 15340420 | Blue <input type="checkbox"/> | | | | | |
| 15350420 | Green <input type="checkbox"/> | | | | | |
| 15370420 | Magenta <input type="checkbox"/> | | | | | |
| 15380420 | White <input type="checkbox"/> | | | | | |
| 15300430 | Blue <input type="checkbox"/> | 10 to 18 | 16 | 25 | 12 x 12.2 x 9.2 | 1 |
| 15310430 | Green <input type="checkbox"/> | | | | | |
| 15330430 | White <input type="checkbox"/> | | | | | |
| 12317629 | Blue <input type="checkbox"/> | | | | | |
| 15340430 | Green <input type="checkbox"/> | | | | | |
| 15360430 | Magenta <input type="checkbox"/> | 15 to 50 | 9 | 30 | 11 x 11 x 8.5 | 1 |

Full-Size Polypropylene


| Cat. No. | Colour | Tubes mL | Wells | Well Dia. mm | L x W x H cm | Qty |
|----------|----------------------------------|----------|-------|--------------|-------------------|-----|
| 15370430 | White <input type="checkbox"/> | 5 to 10 | 72 | 13 | 20.2 x 10.4 x 5.9 | 1 |
| 15380430 | Blue <input type="checkbox"/> | | | | | |
| 15390430 | Green <input type="checkbox"/> | | | | | |
| 15300440 | Magenta <input type="checkbox"/> | | | | | |
| 15310440 | White <input type="checkbox"/> | | | | | |
| 15320440 | Blue <input type="checkbox"/> | 5 to 10 | 72 | 16 | 25 x 12.7 x 7 | 1 |
| 15340440 | Green <input type="checkbox"/> | | | | | |
| 15360440 | Magenta <input type="checkbox"/> | | | | | |

Fisherbrand™ Coated Wire Racks

Chemically resistant HDPE Coated Wire Racks

- Holds all brands and types of conical, round or flat bottom, glass or plastic, tubes and/or bottles
- Square wells and three deck levels help keep tubes secure

Autoclavable

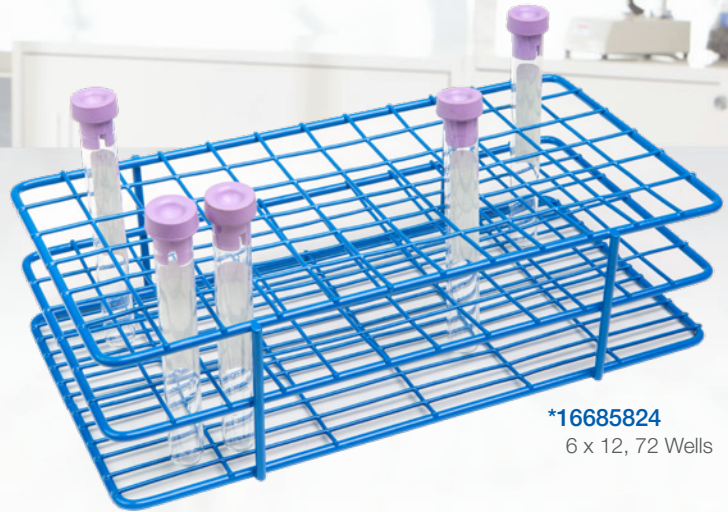

***16685824**
6 x 12, 72 Wells

Fits 10 to 13 mm Tubes

| Cat. No. | Colour | Tube Array | Wells | Well Dia. mm | L x W x H cm | Qty |
|------------------|-----------------|------------|-------|--------------|-------------------|-----|
| 16685824* | Blue | 6 x 12 | 72 | 13 | 20.5 x 10.5 x 6.5 | 1 |
| 16695824 | NEW Blue | 9 x 12 | 108 | 13 | 20.4 x 16.1 x 6.5 | 1 |



***16615834**
4 x 12, 48 Wells

Fits 13 to 16 mm Tubes


| Cat. No. | Colour | Tube Array | Wells | Well Dia. mm | L x W x H cm | Qty |
|------------------|-----------------|------------|-------|--------------|-------------------|-----|
| 16605834 | NEW Blue | 4 x 10 | 40 | 16 | 20 x 8.3 x 5.9 | 1 |
| 16615834 | NEW Blue | 4 x 12 | 48 | 16 | 23.7 x 12.7 x 6.5 | 1 |
| 16625834* | NEW Blue | 9 x 12 | 108 | 16 | 24.1 x 18.8 x 6.5 | 1 |

Fits 16 to 18 mm Tubes – 15 mL

| Cat. No. | Colour | Tube Array | Wells | Well Dia. mm | L x W x H cm | Qty |
|------------------|--------|------------|-------|--------------|------------------|-----|
| 16625834* | Blue | 3 x 5 | 15 | 18 | 12.7 x 7.6 x 7.6 | 1 |

Fits 16 to 20 mm Tubes – 15 mL


| Cat. No. | Colour | Tube Array | Wells | Well Dia. mm | L x W x H cm | Qty |
|-----------------|-----------------|------------|-------|--------------|--------------|-----|
| 16635834 | NEW Blue | 4 x 10 | 40 | 20 | 24 x 9 x 8 | 1 |
| 16645834 | NEW Blue | 8 x 10 | 80 | 20 | 24 x 19 x 8 | 1 |

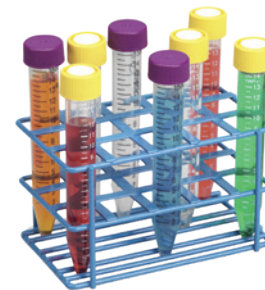

***11749128**
3 x 5, 15 Wells

Fisherbrand™ Coated Wire Racks

Chemically resistant HDPE Coated Wire Racks

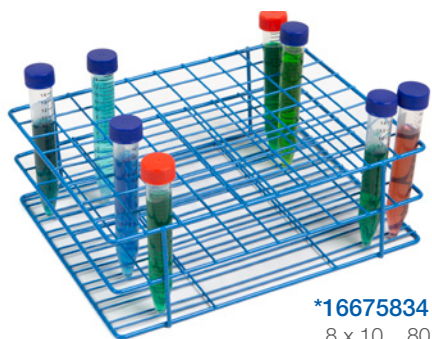

*16675834
8 x 10, 80 Wells

Fits 20 to 24 mm Tubes

| Cat. No. | Colour | Tube Array | Wells | Well Dia. mm | L x W x H cm | Qty |
|----------|-----------------|------------|-------|--------------|---------------|-----|
| 16655834 | NEW Blue | 4 x 10 | 40 | 24 | 26 x 10.7 x 8 | 1 |

Fits 20-25 mm Tubes

| Cat. No. | Colour | Tube Array | Wells | Well Dia. mm | L x W x H cm | Qty |
|-----------|-----------------|------------|-------|--------------|-----------------|-----|
| 16685834 | NEW Blue | 3 x 8 | 24 | 25 | 27.3 x 10.5 x 8 | 1 |
| 16665834 | Blue | 4 x 10 | 40 | 25 | 29 x 11.8 x 8 | 1 |
| 16675834* | Blue | 8 x 10 | 80 | 25 | 29 x 23 x 8 | 1 |

Fits 25 to 30 mm Tubes - 50 mL

| Cat. No. | Colour | Tube Array | Wells | Well Dia. mm | L x W x H cm | Qty |
|-----------|-----------------|------------|-------|--------------|-------------------|-----|
| 11707285* | Blue | 2 x 4 | 8 | 30 | 14.7 x 8 x 8.9 | 1 |
| 16695834 | NEW Blue | 6 x 8 | 48 | 30 | 27.7 x 20.9 x 9.8 | 1 |



*11707285
2 x 4, 8 Wells



*16605844
3 x 8, 24 Wells

Fits 30 to 40 mm Tubes

| Cat. No. | Colour | Tube Array | Wells | Well Dia. mm | L x W x H cm | Qty |
|-----------|-----------------|------------|-------|--------------|-------------------|-----|
| 16605844* | NEW Blue | 3 x 8 | 24 | 40 | 35.7 x 13.9 x 9.8 | 1 |
| 16615844 | NEW Blue | 6 x 8 | 48 | 40 | 35 x 26 x 9.8 | 1 |

Fits 55 to 60 mm Bottles - Centrifuge Bottles

| Cat. No. | Colour | Tube Array | Wells | Well Dia. mm | L x W x H cm | Qty |
|-----------|--------|------------|-------|--------------|--------------------|-----|
| 11805870* | Blue | 2 x 3 | 6 | 60 | 20.5 x 13.5 x 10.4 | 1 |



*11805870
2 x 3, 6 Wells

Fisherbrand™ Reversible 96-Well Tube Racks

Easy-to-grasp tabs for sample handling and storage

Holds 96 x 0.5 mL tubes on one side and 96 x 1.5/2.0 mL on the other.

Optimum Function & Identification

- Wells arranged in standard 8 x 12 array
- Moulded-in grid lines separate 0.5 mL wells making them easier to see
- Raised rings around the 1.5mL/2.0 mL wells make tubes easier to remove
- Alphanumeric identification for quick location and retrieval of tubes

Secure lid

- Translucent lid snaps onto rack tabs for secure, visible storage
- Spacious enough to fit most screw-top micro-centrifuge tubes
- Stackable, saving valuable space

| Cat. No. | Colour | Cat. No. | Colour |
|----------|----------------------------------|----------|--------|
| 11728174 | Assorted* | 11968014 | Green |
| 11738174 | Natural <input type="checkbox"/> | 11978014 | Yellow |
| 11958014 | Blue | | |

*Assorted colours (Blue, Green, Yellow, Orange, Pink)

| Tubes mL | Wells | Well Dia. mm | L x W x H cm | Qty |
|----------|-------|--------------|-----------------|-----|
| 1.5/2 | 96 | 11 | 24.4 x 12 x 5.1 | 5 |
| 0.5 | | 7.6 | | |



Autoclavable

Fisherbrand™ Microtube 80-Well Racks

Popular fraction-collector-plate style rack


Compact 5 x16 array makes filling and transfer easy

- Frosted labelling area on bottom of rack
- Embossed index numbers for quick sample ID and retrieval
- Made from polypropylene

| Cat. No. | Colour |
|----------|----------------------------------|
| 11720344 | Natural <input type="checkbox"/> |
| 11968074 | Blue |
| 11978074 | Green |
| 11988074 | Pink |
| 11998074 | Yellow |
| 11908084 | Orange |
| 11728084 | Assorted* |
| 11918084 | Red |

*Assorted colours (Blue, Green, Yellow, Orange, Pink)

| Tubes mL | Wells | Well Dia. mm | L x W x H cm | Qty |
|----------|-------|--------------|------------------|-----|
| 1.5/2.0 | 80 | 11 | 22.5 x 6.7 x 2.8 | 5 |

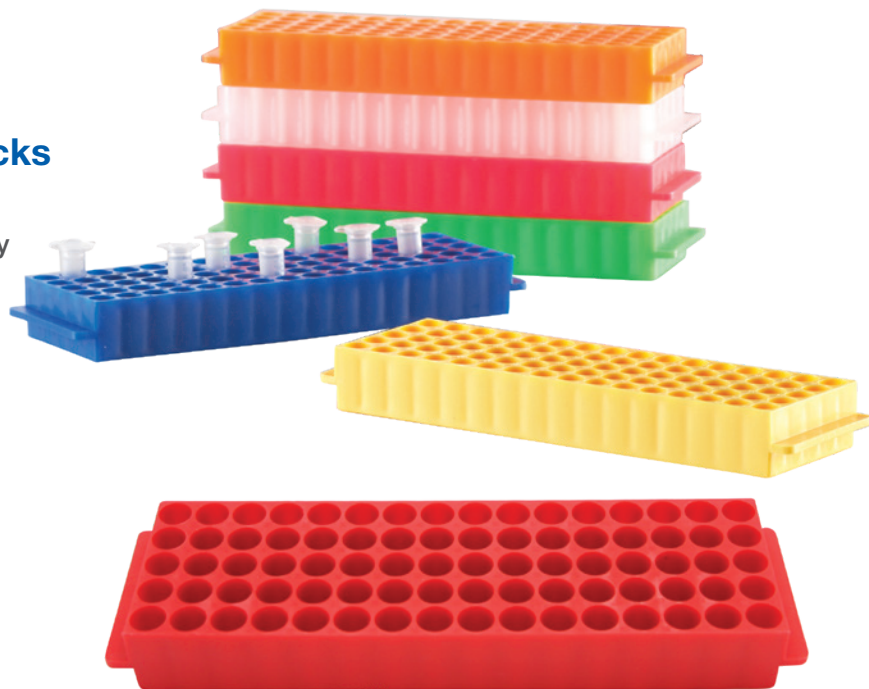


Autoclavable

Fisherbrand™ Microtube Racks

Handles for portability

Durable and simple to use.

- Made from reinforced nylon offering superior chemical and heat-resistance over most other plastic based racks
- Handles provide easy and safe transportation
- Stable design can be used in water baths
- Alphanumeric identification for quick location and retrieval of tubes

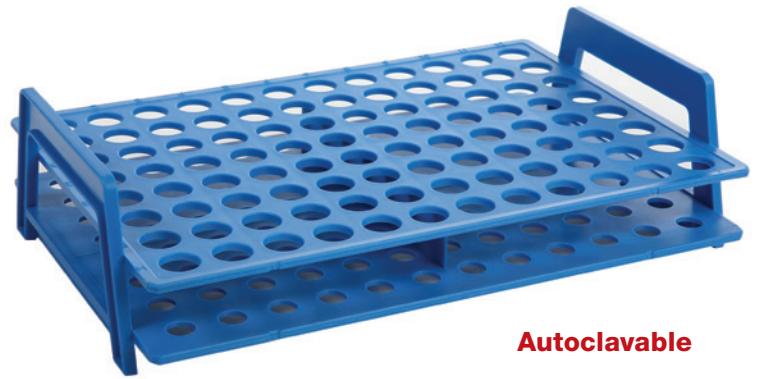

Autoclavable

| Cat. No. | Colour | Tubes mL | Wells | Well Dia. mm | L x W x H cm | Qty |
|----------|--------|----------|-------|--------------|-----------------|-----|
| 15377997 | Blue | 1.5/2.0 | 96 | 11.1 | 21 x 13.1 x 5.4 | 4 |


NEW SIZE

Fisherbrand™ Tube Rack™

High-capacity storage for your tubes

Fully-assembled, polypropylene racks are great for labs requiring large sample storage

- Sturdy, material for easy and secure transport of glass tubes
- Alphanumeric identification for quick location and retrieval of tubes
- Eliminates the need for multiple racks
- Not recommended for autoclaving

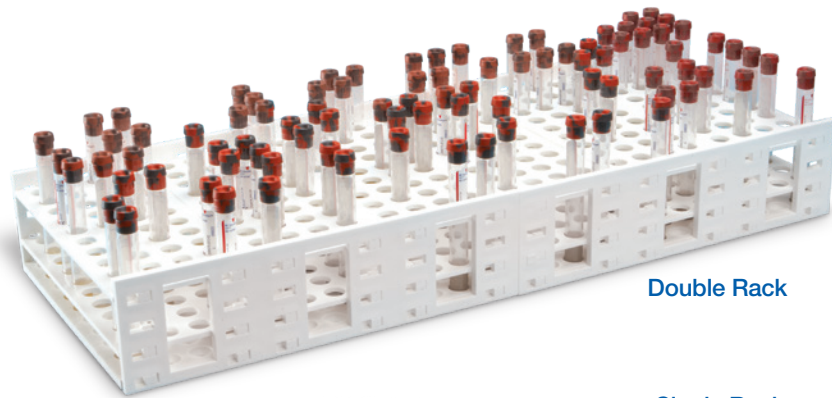


Double Rack

Single Rack

| Cat. No. | Colour | Tubes mL | Wells | Well Dia. mm | L x W x H cm | Qty |
|----------|--------|----------|-------|--------------|-------------------|-----|
| 11720634 | Blue | 5 to 7 | 216 | 13 | 33.7 x 23.5 x 9.1 | 1 |
| 14801442 | White | 10 | 120 | 16 | 33.7 x 26.7 x 9.1 | 1 |
| 11730634 | Blue | 18 to 20 | 120 | 20 | 37.7 x 30.1 x 9.1 | 1 |
| 14811442 | White | | | | | |
| 17157825 | Blue | | | | | |
| 17147825 | White | | | | | |

Double Rack

| Cat. No. | Colour | Tubes mL | Wells | Well Dia. mm | L x W x H cm | Qty |
|----------|--------|----------|-------|--------------|-----------------|-----|
| 15509811 | White | 5 to 7 | 432 | 13 | 67.3 x 23.5 x 9 | 1 |
| 15519811 | Blue | 10 | 240 | 16 | 67.3 x 26.7 x 9 | 1 |
| 15529811 | White | | | | | |
| 15539811 | Blue | | | | | |

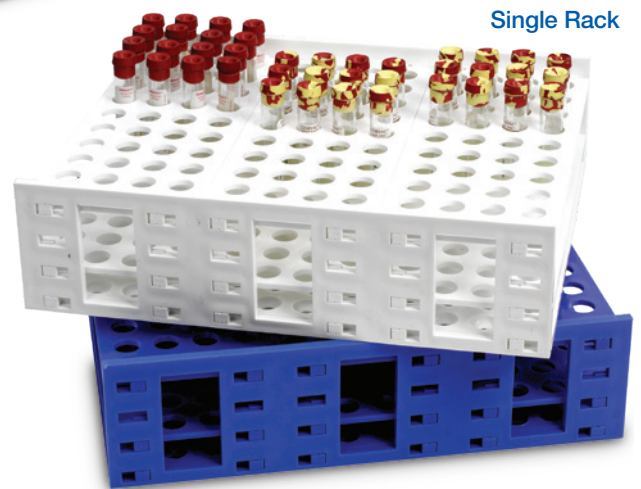

Single Rack

Fisherbrand™ One-Well Friction-Fit Tube Rack

- Sturdy polystyrene base
- Open design enables easy sample identification
- Holds conical- and round- bottom tubes

| Cat. No. | Colour | Tubes mL | Wells | Well Dia. mm | L x W cm | Qty |
|----------|--------|----------|-------|--------------|----------|-----|
| 11354055 | Blue | 50 | 1 | 29 | 9.2 x 3 | 5 |



Fisherbrand™ Expanding Storage Box

Single rack, double function

This convenient working and storage rack can help reduce sample handling time and eliminates the need for multiple racks when testing samples kept at cold temperatures.

Efficient Handling

- Expand rack instantly creating more working space without interference from open microfuge lids
- Collapse rack into closed position for freezer storage
- Compact design fits into standard 53mm freezer racking system
- Alphanumeric identification for quick location and retrieval of tubes

Safe and Stackable Storage

- Clear lids locks the rack together for safe handling and transportation
- Low profile, easily stackable

Expanding Storage Box

- Holds 64 x 1.5/2.0 mL cryo and/or microtubes
- Expanding wells have a centre depression to stabilise conical tubes and a cryogenic vial lock for one handed cap removal



Durable Construction

Sides made from polycarbonate, Tube holders made of polypropylene

Autoclavable



Clear Lid

Locks rack together for safe handling and transportation



Expanding Storage Box

| Cat. No. | Colour | Tubes mL | Wells | Well Dia. mm | L x W x H (Open) cm | L x W x H (Closed w/lid) cm | Qty |
|----------|--------|----------|-------|--------------|---------------------|-----------------------------|-----|
| 15162855 | White | 1.5/2.0 | 64 | 12.3 | 31 x 12.8 x 3.4 | 13.3 x 13.3 x 5.3 | 1 |



Fisherbrand™ Tube Rack

Convenient Rack System

This storage box with eight removable tube racks lets you easily remove an individual sample or an entire segment and leave the rest.

- Printed index letters on each segment make for quick sample ID and retrieval
- Moisture and mould resistant
- Thin film polypropylene boxes are 51 mm high and fit inside standard freezer racks
- Boxes and racks can be stored down to temperatures as low as -80°C
- Eight polypropylene segments hold 1.5–2.0 mL cryo and microcentrifuge tubes (64 tubes)
- Wells feature a centre depression to stabilise conical tubes and a cryogenic vial lock for one-handed cap removal



Autoclavable

Tube Rack

| Cat. No. | Colour | Colour | Tubes mL | Wells | Well Dia. mm | L x W x H cm | Qty |
|----------|--------|--------|----------|-------|--------------|-------------------|-----|
| 15346548 | Purple | ■ | 1.5/2.0 | 64 | 12.3 | 12.9 x 12.9 x 3.5 | 1 |

Microtube Storage Boxes

Stylish & Identifiable

- Made of translucent, vibrantly coloured polypropylene for visibility and colour coding
- Tube capacity is clearly indicated on the outside of the box
- Uniform wells and alphanumeric location identification and a frosted labelling area

Robust

- Stackable, saving valuable space in the lab
- Can withstand temperatures as low as -80°C

Autoclavable



Fisherbrand™ 50-Well Microtube Storage Boxes

Ideal for sample or case study storage

- For smaller storage requirements of 1.5/2.0 mL microtubes or cryogenic vials.
- Easy-open hinged lid with snap clasp

| Cat. No. | Colour |
|----------|----------------------------------------------------------------------------------------------------------------------------------------|
| 11704316 | Natural <input type="checkbox"/> |
| 10707992 | Assorted* <input type="checkbox"/> <input type="checkbox"/> <input type="checkbox"/> <input type="checkbox"/> <input type="checkbox"/> |

*Assorted colours (Blue, Green, Purple, Yellow, Orange)

| Tubes mL | Wells | Well Dia. mm | L x W x H cm | Qty |
|----------|-------|--------------|------------------|-----|
| 1.5/2.0 | 50 | 12.3 | 14.2 x 9.2 x 5.6 | 5 |

Fisherbrand™ 81-Well Microtube Storage Boxes

Compact footprint for microtube and cryogenic vial storage

- Store 1.5/2.0 mL microtubes or cryogenic vials
- Polypropylene boxes with easy open, friction-fit lids
- Black, light blocking version available for light sensitive samples

| Cat. No. | Colour | Cat. No. | Colour |
|----------|----------------------------------|----------|----------------------------------------------------------------------------------------------------------------------------------------|
| 10500203 | Natural <input type="checkbox"/> | 10325002 | Yellow <input type="checkbox"/> |
| 10273222 | Blue <input type="checkbox"/> | 10385042 | Orange <input type="checkbox"/> |
| 10050182 | Green <input type="checkbox"/> | 15558319 | Black <input type="checkbox"/> |
| 10070182 | Pink <input type="checkbox"/> | 10243272 | Assorted* <input type="checkbox"/> <input type="checkbox"/> <input type="checkbox"/> <input type="checkbox"/> <input type="checkbox"/> |

*Assorted colours (Blue, Green, Yellow, Orange, Pink)

| Tubes mL | Wells | Well Dia. mm | L x W x H cm | Qty |
|----------|-------|--------------|---------------|-----|
| 1.5/2.0 | 81 | 12.3 | 13 x 13 x 4.7 | 5 |



Fisherbrand™ 100-Well Microtube Storage Boxes

Secure storage for microtube and cryogenic vials

- For larger storage requirements of 1.5/2.0 mL microtubes or cryogenic vials
- Black, light blocking version available for light sensitive samples
- Easy-open hinged lid with snap clasp

| Cat. No. | Colour | Cat. No. | Colour |
|----------|----------------------------------------------------------------------------------------------------------------------------------------|----------|---------------------------------|
| 11790334 | Natural <input type="checkbox"/> | 15579811 | Green <input type="checkbox"/> |
| 11918014 | Black <input type="checkbox"/> | 15589811 | Yellow <input type="checkbox"/> |
| 11700344 | Assorted* <input type="checkbox"/> <input type="checkbox"/> <input type="checkbox"/> <input type="checkbox"/> <input type="checkbox"/> | 15459123 | Orange <input type="checkbox"/> |
| 15569811 | Blue <input type="checkbox"/> | | |

*Assorted colours (Blue, Green, Yellow, Orange, Pink)

| Tubes mL | Wells | Well Dia. mm | L x W x H cm | Qty |
|----------|-------|--------------|-------------------|-----|
| 1.5/2.0 | 100 | 12.3 | 14.1 x 15.1 x 5.7 | 5 |



Fisherbrand™ Low-Temp PCR Rack: 96-Well

Slide, set, work & store PCR® samples

Multi-purpose rack converts from compact workstation to lidded storage solution accommodating individual and strip 0.2mL tubes

Versatile & Robust Solution

- Polycarbonate lid serves as cover and elevated rack base - just slide open and pivot lid under rack for compact workstation
- Detachable lid can be used in right or left handed orientation
- Stackable saving valuable bench-top and freezer space
- Multiple colours offered for coding options
- Store down to -80°C with appropriate temperature viable tubes

US Patent Pending, Community Design No: 004101756-0001-0002



| Cat. No. | Colour | Cat. No. | Colour |
|----------|-----------|----------|--------|
| 15758177 | Assorted* | 15778177 | Green |
| 15768177 | Blue | 15788177 | Purple |

*Assorted colours (Purple, Blue, Green, Orange, Yellow)

| Tubes mL | Wells | Well Dia. mm | L x W x H cm | Qty |
|----------|-------|--------------|------------------|-----|
| 0.2 | 96 | 5.2 | 12.5 x 9.2 x 3.7 | 5 |

Fisherbrand™ PCR Tube Racks 96-Well

Perfect for pre- or post-PCR® sampling

- Made from polypropylene
- Accommodates individual tubes or strips
- Removable hinged lid lies flat
- Alphanumeric identification for quick location and retrieval of tubes

Autoclavable

| Cat. No. | Colour | Cat. No. | Colour |
|----------|-----------|----------|--------|
| 11710344 | Assorted* | 11948074 | Blue |
| 11527593 | Natural | 11958074 | Green |

*Assorted colours (Blue, Green, Pink, Yellow, Orange)

| Tubes mL | Wells | Well Dia. mm | L x W x H cm | Qty |
|----------|-------|--------------|------------------|-----|
| 0.2 | 96 | 5.2 | 12.7 x 9.9 x 3.0 | 5 |



Fisherbrand™ PCR Tube Racks Reversible

Removable hinged lid fits both sides

- Made from polypropylene
- PCR side of the rack holds 8 or 12 x 0.2 mL tube strips or individual tubes
- Reverse side holds 0.5 mL tubes
- Both sides of the rack have 12 wells for 1.5 mL tubes

Autoclavable

| Cat. No. | Colour | Tubes mL | Wells | Well Dia. mm | L x W x H cm | Qty |
|----------|--------|----------|-------|--------------|-------------------|-----|
| 11787265 | | 0.2 | 168 | 5.5 | 20.6 x 13.1 x 5.3 | 5 |
| | | 0.5 | 40 | 7.6 | | |
| | | 1.5/2.0 | 12 | 11.4 | | |

*Assorted colours (Blue, Green, Pink, Yellow, Orange)



Fisherbrand™ PCR Cooler, 96-Well


NEW

Ice-Free Cooling System For PCR Preparation and Storage
 Cools industry standard format PCR tubes, strips and single format as well as most 96-Well plates, in both standard and low profile well shape

Consistent & Reliable

- Modern design cools faster in freezers and stays cold longer on the bench top
- Maintains the entire 96-well plate at a temperature of 0°C for over an hour
- Maintains sample temperature of less than 5°C for up to three hours (with lid off)
- When the temperature ceiling of 6°C has been exceeded, users are signaled by a colour change of the cooler or a single well

96-Well 8 x 12 configuration

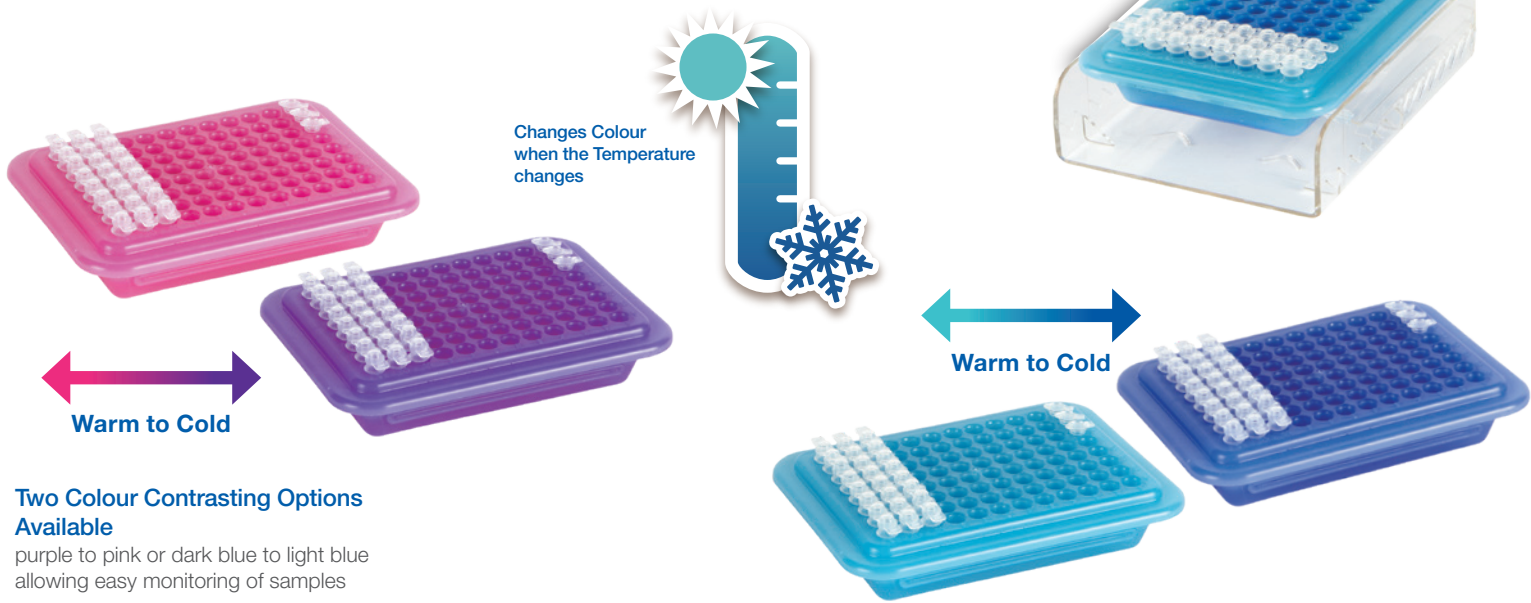

Two Colour Contrasting Options Available

purple to pink or dark blue to light blue allowing easy monitoring of samples



Ergonomic Textured Finger Grips

provide a secure hold, reducing the chance of dropping the PCR Cooler, and makes handling and transportation easy, especially when using gloves



Vented Lid

enables faster cooling of samples while in freezer



Lid doubles as a stand

When used as a stand the space between the cooler and bench top enables the cooler to stay cold longer



Alphanumeric Well Identification

for quick location and retrieval of tubes

| Cat. No. | Warm Colour | Cold Colour | Wells | Tube Sizes mm | L x W x H cm | Qty |
|----------|-----------------------------------------------------|---------------------------------------------------|-------|-----------------------|-----------------|-----|
| 16300412 | Pink ■ | Purple ■ | 96 | 0.2, 0.5 Assay Plates | 11 x 14.8 x 4.8 | 2 |
| 16310412 | Light Blue ■ | Dark Blue ■ | | | | |
| 16320412 | Assorted* | | | | | |

*Assorted colours (Blue, Pink)

32 Place PCR Rack

Designed for pre- and post-PCR sampling, streamline your PCR workflow with the 32-Place PCR Rack. This versatile rack accommodates 0.2 mL tubes both individually and in strips of 8, as well as 32-place assay plates (8x4)

Space Saving Design

- Small and compact design is stackable, saving valuable space in the freezer
- Right sized when processing smaller batches
- Lid hinge design allows lid to sit at 180 degrees, clear of sample processing, or can be removed entirely

Reusable Convenience

- Transparent lid protects samples and allows for quick assessment of contents
- Embossed numbers and letters on the wells provide easy identification and tracking of samples
- Light-weighted construction made durable for reusability
- **Autoclavable** at 121° for 15 minutes

| Cat. No. | Tubes mL | No. of Wells | L x W x H in | L x W x H cm | Qty |
|----------|-------------|--------------|-----------------|------------------|-----|
| 17862455 | 0.2 | 32 | 3.8 x 6.9 x 1.2 | 9.8 x 17.5 x 3.1 | 10 |

Stores down to -80°C



Fisherbrand™ Cool Containers

Durable containers for chilling temperature-sensitive samples. Heavy insulation properties for the ultra cold, cool containers are ideal for ice, ice-salt slurries and more

- Made from single-moulded polyurethane
- Close-fit lid keeps evaporation to a minimum and features a comfortable recessed knob
- Lid nests under bucket
- For low-temperature use to -196°C at 1 bar
- Stackable
- Pour spout and large handles

| Cat. No. | Style | Capacity | L x W x H mm | Colour | Qty |
|----------|--------|----------|-------------------|--------|-----|
| 11314085 | Bucket | 4 L | 181 (H) x 234 (D) | Black | 1 |
| 11324085 | Pan | 1 L | 267 x 195.5 x 114 | Green | |
| 11334085 | Pan | 4 L | 381 x 254 x 140 | Blue | |
| 11344085 | Pan | 9 L | 521 x 292 x 163 | Red | |



Lid fits below
bucket to save
space



Fisherbrand™ Foam Tube Racks

High storage, minimal footprint

Two-piece racks ideal for long-term studies or tube storage

- Made from lightweight expanded polystyrene (EPS) foam
- Stack securely with nesting features on lid and base
- Extra spacing between wells allows for easy gripping with fingertips
- Store down to temperatures of -80°C, when used with temperature viable tubes


| Cat. No. | Colour | Tubes mL | Wells | Well Dia mm | L x W x H cm | Qty |
|----------|--------|-------------|-------|----------------|------------------|-----|
| 11774016 | White | 1.5/2.0 | 100 | 11.5 | 33.6 x 9.5 x 7.3 | 10 |
| 11948084 | White | 1.5/2.0 | 50 | 11.5 | 21 x 11 x 7.1 | 10 |



Fisherbrand™ Floating Microtube Rack

For use on bench and in ice and water baths

- Made from polypropylene
- Attach legs for upright use; detach legs to force samples up for easier removal

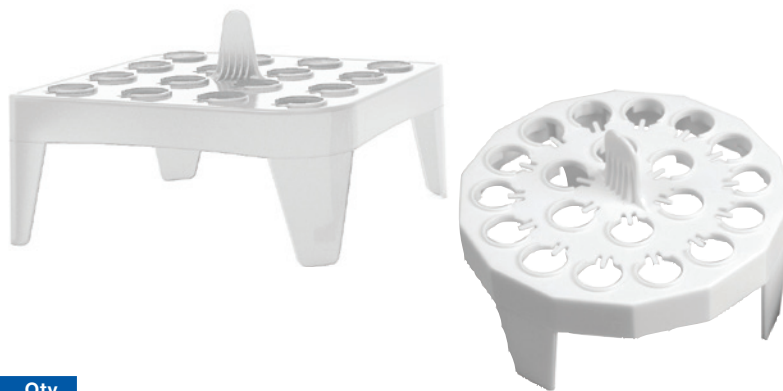

| Cat. No. | Colour | Tubes mL | Wells | Well Dia. mm | L x W x H cm | Qty |
|----------|--------------------------------|----------|-------|--------------|--------------|-----|
| 13436959 | White <input type="checkbox"/> | 1.5/2.0 | 16 | 10.9 | 10.2 | 4 |
| 13456959 | | | 20 | | 9.8 dia. | |

Fisherbrand™ Floating Tube Racks

Reusable and flexible for a variety of applications

- Ideal for immersing tubes in water baths, tubs or beakers
- Brightly-coloured HDPE foam racks resist water absorption
- Racks include detachable carrying handles—easily lift samples out of baths



| Cat. No. | Colour | Tubes mL | Wells | Well Dia. mm | L x W x H cm | Qty |
|----------|--------------------------------------------|-----------------|-------|--------------|-------------------|-----|
| 11704736 | Green <input checked="" type="checkbox"/> | 50 | 4 | 29 | 16.5 x 13.7 x 6.5 | 5 |
| 11714736 | Blue <input checked="" type="checkbox"/> | 15 | 8 | 16 | 17.2 x 5.7 x 6.5 | |
| 11724736 | Yellow <input checked="" type="checkbox"/> | 1.5/2.0 | 24 | 9.3 | 15 x 8.3 x 7 | |
| 11736565 | Blue <input checked="" type="checkbox"/> | 0.2/0.5/1.5/2.0 | 18 | 9.3 | 8.2 dia. | |

"X" Design
to hold a variety of tube sizes

NEW



Sliding Lock
provides secure storage with easily identified lock and unlock symbols.

Slide Rail Design
ensures a secure fit allowing easy placement and removal of slides

Centrally Positioned Numbering
enables simple slide organisation with easy viewing of location numbers

Patent Pending

Fisherbrand™ 100-Place Premium Plus Slide Box

Modern design for maximum functionality

Optimal Organisation

- Organise by "project" and "date" moulded into lid
- Inventory card index corresponds to slide numbers and can be removed to let you write on a flat surface
- Vibrant translucent colours are easy to identify and enable coding options

Functional Design

- Store vertically or horizontally to efficiently use space
- Moulded feet allow horizontal stacking of multiple units
- Unique side storage accommodates a pen, desiccant or other small items
- Hinged lid smoothly opens 180° to lay flat on work surfaces

Safe & Secure

- EVA foam lining helps hold slides securely in place and protect against damage
- Stylish wave surface texture is easy to grip for safe transport

100-Place Premium Plus Slide Box

| Cat. No. | Colour | Cat. No. | Colour |
|----------|--------|----------|-----------|
| 15810086 | Blue | 15840086 | Red |
| 15820086 | Orange | 15850086 | Purple |
| 15830086 | Green | 15860086 | Assorted* |

*Assorted colours (Purple, Blue, Green, Red, Orange)

| Lining | Hinge | L x W x H cm | Qty |
|----------|---------------|-------------------|-----|
| EVA Foam | Sliding Latch | 20.7 x 19.5 x 3.2 | 5 |

Inventory Card

| Cat. No. | L x W x H cm | Qty |
|----------|--------------|-----|
| 15870086 | 21.6 x 28 | 10 |

Fisherbrand™ Premium Microscope Slide Boxes

Slide boxes are perfect for storage or transportation of standard size slides 26 x 76 mm (1 x 3 inch)

Organise & Identify

- Lining and slots are numbered to correspond to the slide inventory sheet
- Organise by "Project" and "Date" with pre-moulded labelling area
- Available in multiple colours for project coding

Durable & Safe

- ABS plastic boxes available with either cork or foam lining
- Designed to securely store slides in grooved slots
- Moulded feet secure stacking for easy transport and storage surfaces and prevents misplacement



25-Place Slide Box

Easily access slides with friction-fit cover



| Cat. No. | Colour | Lining | Hinge | L x W x H cm | Qty |
|----------|--------------------------------|--------|--------------|------------------|-----|
| 11354135 | White <input type="checkbox"/> | Cork | Friction Fit | 14.1 x 8.8 x 3.5 | 1 |
| 11314135 | Blue <input type="checkbox"/> | | | | |
| 11324135 | Green <input type="checkbox"/> | | | | |
| 11334135 | Red <input type="checkbox"/> | | | | |
| 11711476 | Blue <input type="checkbox"/> | Foam | Friction Fit | 14.1 x 8.8 x 3.5 | 1 |
| 11324135 | Green <input type="checkbox"/> | | | | |
| 11721476 | Red <input type="checkbox"/> | | | | |
| 11731476 | White <input type="checkbox"/> | | | | |

50-Place Slide Box

Thumb-latch lock mechanism



| Cat. No. | Colour | Lining | Hinge | L x W x H cm | Qty |
|----------|--------------------------------|--------|-------------|------------------|-----|
| 11731486 | Blue <input type="checkbox"/> | Cork | Thumb-Latch | 20.9 x 8.6 x 3.5 | 1 |
| 11771486 | Red <input type="checkbox"/> | | | | |
| 11781486 | White <input type="checkbox"/> | | | | |

100-Place Slide Box

Closes securely with a rust-resistant nickel-plated clasp and pin hinges



| Cat. No. | Colour | Lining | Hinge | L x W x H cm | Qty |
|----------|---------------------------------|--------|--------------|-------------------|-----|
| 11314145 | Blue <input type="checkbox"/> | Cork | Nickel Clasp | 20.8 x 17.5 x 3.4 | 1 |
| 11374145 | Black <input type="checkbox"/> | | | | |
| 11324145 | Green <input type="checkbox"/> | | | | |
| 11334145 | Red <input type="checkbox"/> | | | | |
| 11344145 | Yellow <input type="checkbox"/> | | | | |
| 11354145 | White <input type="checkbox"/> | Foam | Nickel Clasp | 20.8 x 17.5 x 3.4 | 1 |
| 11364145 | Gray <input type="checkbox"/> | | | | |
| 11701486 | Blue <input type="checkbox"/> | | | | |
| 10428241 | Green <input type="checkbox"/> | | | | |
| 11711486 | Red <input type="checkbox"/> | | | | |
| 11394135 | Yellow <input type="checkbox"/> | Cork | Nickel Clasp | 20.8 x 17.5 x 6 | 1 |
| 11721486 | White <input type="checkbox"/> | | | | |

100-Place Tall Slide Box

Additional space for larger 75 mm x 50 mm. slides
Stackable with 100-Place Premium Slide Boxes



| Cat. No. | Colour | Lining | Hinge | L x W x H cm | Qty |
|----------|-------------------------------|--------|--------------|-----------------|-----|
| 15186823 | Blue <input type="checkbox"/> | Cork | Nickel Clasp | 20.8 x 17.5 x 6 | 1 |

Fisherbrand™ Standard Microscope Slide Boxes

Constructed of durable polypropylene

Improved Identification

- Moulded slot numbers correspond directly to the slide inventory card
- Organise with "Project" and "Date" pre-moulded labelling area

Secure & Stackable

- Hinged lids are held securely closed with a latch lock
- Moulded feet secure stacking for easy transport and storage

25-Place Slide Box

| Cat. No. | Colour | Lining | Hinge | L x W x H cm | Qty |
|----------|--------|-------------------------------------|-----------------|----------------|-----|
| 11741476 | Blue | <input checked="" type="checkbox"/> | None Latch Lock | 14.1 x 9.2 x 3 | 1 |
| 11751476 | Red | <input checked="" type="checkbox"/> | | | |
| 11761476 | White | <input type="checkbox"/> | | | |

100-Place Slide Box

| Cat. No. | Colour | Lining | Hinge | L x W x H cm | Qty |
|----------|--------|-------------------------------------|-----------------|-----------------|-----|
| 11771476 | Blue | <input checked="" type="checkbox"/> | None Latch Lock | 21 x 16.9 x 3.7 | 1 |
| 11781476 | Red | <input checked="" type="checkbox"/> | | | |
| 11791476 | White | <input type="checkbox"/> | | | |



Holds standard slides
26 x 76 mm



Thumb Notch
allows for easy removal of slides

Fisherbrand™ Cardboard Slide Tray

Cardboard slide trays hold standard slides
26 x 76 mm

Cardboard slide folders offer a convenient solution for safely working, transporting and short-term storing of slides

Organise & Protect

- Index on the front cover corresponds with the number grid inside for identification
- Recessed slide compartments and inside covers for specimen protection
- Slides cannot be accessed when cover is closed reducing contamination

Sturdy Construction

- Compartment holder is made from a double layer of densely pressed cardboard. Cover is single layer cardboard
- Front protection flaps are attached with durable tape

| Cat. No. | Colour | Capacity | L x W x H cm | Qty |
|----------|--------|-----------|-------------------|-----|
| 11703217 | Black | 20 Slides | 20.8 x 34.1 x 0.9 | 1 |
| 11790465 | Blue | | | |
| 12942017 | Green | 30 Slides | 28.1 x 32.8 x 1 | 1 |
| 11713217 | Black | | | |

NEW

Fisherbrand™ Slide Folder Clasps

Keep slide folder flaps securely closed

- Ideal for histology, microbiology and pathology as well as multiple user research environments
- Clasps can be labeled with a permanent marker or affixed label




| Cat. No. | Colour | L x W x H cm | Qty |
|----------|----------------------------------------------------------------------------------------------------------------------------------------|-----------------|-----|
| 15813977 | Assorted* <input type="checkbox"/> <input type="checkbox"/> <input type="checkbox"/> <input type="checkbox"/> <input type="checkbox"/> | 5.6 x 4.4 x 1.6 | 100 |
| 15823977 | Clear <input type="checkbox"/> | | |

*Assorted colours (Clear, Gray, Blue, Green, Red)

Fisherbrand™ Slide Holders

Hold slides during evaluation, cleaning and drying

- Made from chemical and stain-resistant polyethylene
- Solid, strong construction with rubber feet for benchtop stability
- Angled slide holder slots allow easy access to individual slides and viewing of slide labels

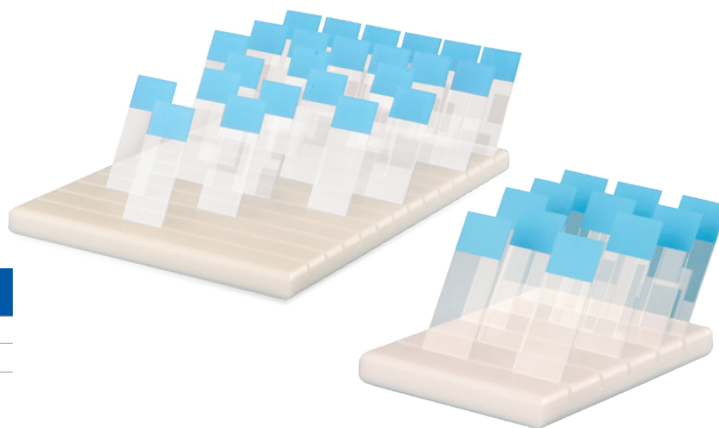

| Cat. No. | Colour | Capacity | L x W x H cm | Qty |
|----------|----------------------------------|-----------|-------------------|-----|
| | Natural <input type="checkbox"/> | 40 Slides | 29.8 x 10.5 x 2.4 | 1 |
| | Natural <input type="checkbox"/> | 78 Slides | 29.9 x 17.8 x 2.1 | 1 |

Fisherbrand™ Slide Staining Tray Set

Reduce possible cross-contamination and messy work areas. Rugged staining tray is an economical alternative to metal versions

- Made from a polypropylene and polyethylene blend for stability
- Stain, rinse and dry slides on a single working tray
- Fits up to eight slides and deep well holds up to 38 mL
- Compact size and recessed handles for easy transport from work area to sink. Pack of 4 base trays and 1 lid



| Cat. No. | Colour | Capacity | L x W x H cm | Qty |
|----------|-------------------------------------------|----------|-----------------|-----|
| 11968004 | Black <input checked="" type="checkbox"/> | 8 Slides | 30.5 x 12.7 x 3 | 4 |

Fisherbrand™ 18-Place Slide Staining Tray

Ideal for clinical & research laboratory staining

Tray holds standard microscope slides (76 x 26 mm) securely in place for staining, rinsing and drying applications

Robust design

- Durably constructed from solvent resistant polypropylene blend
- Non-skid Silicone feet keep tray in place

Convenient to use

- Two ways to empty stain liquid
 - Pour spouts on each corner
 - rain hole with plug
- Equipped with handle cut-outs for easy handling

Versatile for multiple testing formats

- Black lid included for light sensitive applications
- Humidity retention wells for specialized staining needs

NEW



Dark Lid
protects slides for light-sensitive applications

Convenient Pour Spouts
for quick disposal of liquid waste



18-Place Slide Staining Tray

| Cat. No. | Color | Capacity | L x W x H in | L x W x H cm | Qty |
|----------|-------------------------------------------|----------|---------------------|------------------|----------------|
| 17784749 | Black <input checked="" type="checkbox"/> | 18-Place | 16.1 x 8.9 x 1.6 in | 41 x 22.5 x 4 cm | 1 ea, 10 ea/cs |

Fisherbrand™ Slide Mailers

Safe, secure transport of microscope slides

Cardboard Slide Mailers

| Cat. No. | Colour | Capacity | L x W x H mm | Qty |
|----------|---------------------------------------------|----------|---------------|-----|
| 12323107 | Natural <input type="checkbox"/> | 1-Place | 95 x 42 x 5 | 25 |
| 15873967 | Natural <input checked="" type="checkbox"/> | 4-Place | 120 x 140 x 6 | 20 |

- Heavy cardboard mailers for safe transport
- Cover that opens for viewing and thumb grooves for easy slide removal
- Securely holds slides with or without cover glass
- Flat surface area is ideal for labelling or writing necessary information
- Each will hold 26 x 76 mm (1 x 3 in)



Plastic 2-Place Slide Mailers

| Cat. No. | Colour | Capacity | L x W x H mm | Qty |
|----------|-----------------------------------------------|----------|--------------|-----|
| 11745088 | Assorted* <input checked="" type="checkbox"/> | 2-Place | 93 x 68 x 5 | 25 |
| 11755088 | Natural <input type="checkbox"/> | 1-Place | 83 x 41 x 5 | 50 |
| 15863967 | Natural <input type="checkbox"/> | 1-Place | 83 x 41 x 5 | 50 |

*Assorted colours (Blue, Purple, Red, Orange, Yellow)

- Ideal for secure shipping and storing of slides
- Polypropylene protectors hold two standard (26 x 76 mm (1 x 3 in) slides securely
- Large, easy-to-secure snaps provide safe storage for valuable slides
- Variety of colours brings some fun to the lab and enables colour coding



Plastic 5-Place Slide Mailers

| Cat. No. | Colour | Capacity | L x W x H mm | Qty |
|----------|-----------------------------------------------|----------|--------------|-----|
| 11714316 | Natural <input type="checkbox"/> | 5-Place | 81 x 19 x 31 | 25 |
| 11719885 | Natural <input type="checkbox"/> | 5-Place | 43 x 24 x 88 | 25 |
| 16310736 | Assorted* <input checked="" type="checkbox"/> | 5-Place | 43 x 24 x 88 | 25 |

*Assorted colours (Blue, Purple, Red, Orange, Yellow)

- Ideal for secure shipping, storing or staining slides
- Polypropylene containers hold five standard slides 26 x 76 mm (1 x 3 in)
- Inside wall features slotted channels to keep slides separated and to make removal easier
- Lid closes securely to keep slides safe





Complete system for easy counting

Fisherbrand™ Colony Counter

Rapid accurate colony counter

Combines a multi-function electronic counter with a Sharpie™ pen for marking to prevent missing or double-counting colonies

- Light, ergonomic body is balanced for easy manipulation
- Select between an audible beep or an LED flash to acknowledge count
- Store up to 32 separate counts and display the total of all counts
- Count up/down verification
- Quick and easy to change between different coloured markers to check different colony counts in the same view
- Lock counter off to allow use of the marker without registering more counts

| Cat. No. | Description | Qty |
|----------|---------------------------|-----|
| 15253305 | Electronic Colony Counter | 1 |



Weighted stand for bench top use



Thumb button allows you to use as a digital tally counter

Package includes:

- Colony Counter
- Black Sharpie® pen for highly visible marks on both plastic and glass
- Long-life button battery
- Hand tally adapter counts without marking
- Weighted stand for convenient bench top use



Fisherbrand™ Mini-Magnifier

Perfect accessory for the Colony Counter

- Swiveling, moulded lens. 1.75x magnification
- Enlarges objects for easier viewing

| Cat. No. | Description | L x W x H mm | Qty |
|----------|----------------|--------------------|-----|
| 12816355 | Mini Magnifier | 161 x 135 x 207 mm | 1 |

Fisherbrand™ Benchtop Counters

Mechanical, benchtop counters are perfect for blood cells (label included), bacterial colonies or other repetitive counting

- Each key records up to 999
- Total window keeps track of total number of strokes made on all other keys; bell rings when total reaches 100
- Blood cell key labels



| Cat. No. | No. of Keys | L x W x H cm | Qty |
|----------|-------------|----------------|-----|
| 11334095 | 6-Key | 25 x 7.6 x 5.6 | 1 |
| 11344095 | 9-Key | 32 x 7.6 x 5.6 | 1 |

Fisherbrand™ Colony Counter Bundle

Quickly identify and tally colonies with precision using this all-inclusive package designed to streamline all colony counting procedures.

Fisherbrand eCount Colony Counter

Combines an electronic counter with a Sharpie™ pen for marking to prevent missing or double-counting colonies. Features count verification with beeps or LED flashes, dual counting modes (up and down), a display for total counts (up to 32 separate counts), and the ability to disable the counter for writing totals and notes.

Fisherbrand Lightbox

The Fisherbrand Lightbox serves as a LED light source for efficient bacterial colony counting. It features a spacious viewing area (125 x 100 mm) suitable for various samples, and is powered by six 1.5V AAA batteries (included).

Fisherbrand Magnifier

The Fisherbrand Mini-Magnifier enhances object visibility by providing 1.75X magnification. It features a swiveling, molded lens for comfortable observation.

Fisherbrand Silicone Lab Mat

This FDA food-grade silicone labmat is suitable for various laboratory settings, including cold rooms, clean rooms, and sterile suites, where high sterility is essential. Effectively safeguards countertops from temperatures up to 200°C (392°F) and demonstrates general chemical resistance to moderate oxidizing chemicals.

- **Autoclavable** at 121° for 15 minutes



Reversible
Blue and Gray

| Cat. No. | Description | Qty |
|----------|------------------------------------|-----|
| 17842966 | Fisherbrand Colony Counter Bundle* | 1 |

Fisherbrand eCount Colony Counter (15253305)

| L x W x H in | L x W x H cm |
|-----------------|-----------------|
| 7.4 x 1.2* | 17 x 3* |

*Measurements shown are with Sharpie® Pen inserted

Package Includes

| | | |
|-----------------------------------|--------------------|--------------------|
| eCount Colony Counter | Black Sharpie® Pen | Hand tally Adapter |
| Long-life button battery included | Weighted Stand | |





Fisherbrand Lightbox & Magnifier (16360666, 12816355)

| Description | L x W x H in | L x W x H cm |
|-----------------|-----------------|--------------------|
| Mini-Light Box* | 6.2 x 5.9 x 0.6 | 15.8 x 15 x 1.4 |
| Mini-Magnifier | 6.3 x 5.3 x 8.1 | 16.1 x 13.5 x 20.7 |

Fisherbrand Silicone Lab Mat (15837944)

| Colour | L x W x H in | L x W x H cm |
|-----------|-----------------|-----------------|
| Grey/Blue | 14 x 23.5 x 0.1 | 35 x 60 x 0.2 |

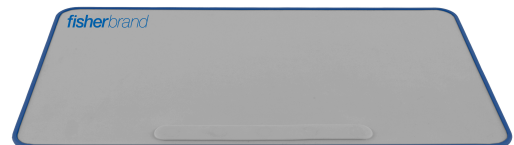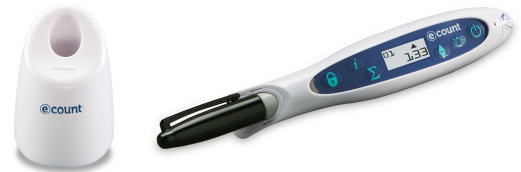



Ratcheting Handle
folds out of the way for easy access,
even over tall sharps containers



Supply Drawer

Side Pockets

Fisherbrand™ Blood Collection Tray

Designed with phlebotomists in mind

- Keeping a clean and organised tray has never been easier. Wash the ABS tray as needed, and replace the insert periodically for a fresh start
- Multiple compartments keep a variety of supplies neatly at hand and organised
- Built in side pockets offer additional storage space
- Handle includes soft rubber grip to reduce hand and finger strain
- Each tray includes either a 13 or 16 mm tube rack and two disposable polystyrene inserts

Configurations

| Cat. No. | Tray Insert | Tube Rack mm | L x W x H cm | Handle Height cm | Qty |
|----------|-------------|--------------|--------------------|------------------|-----|
| 11740634 | 11760634 | 13 | 42.9 x 26.5 x 17.2 | 28 | 1 |
| 11750634 | 11770634 | 16 | | | |

Insert Trays

| Cat. No. | Description | L x W x H cm | Qty |
|----------|---------------------------|--------------|-----|
| 11760634 | Insert A for 13 Tube Rack | 36 x 21 | 4 |
| 11770634 | Insert B for 16 Tube Rack | | |

Tray Inserts
come in two configurations and
include an extra-deep compartment
for a sharps container.



Disposable Insert A
for 13 mm tube rack



Disposable Insert B
for 16 mm tube rack

US Patent No. D544703



Keeps contents clean and secure



Fisherbrand™ Tray Lid

Cover for the droplet Blood Collection Tray

Protect the contents of the Blood Collection Tray with custom designed, secure lid

- Notched design provides access to tray handle to maintain portability
- Made of easy to clean polystyrene
- Can be secured with cable ties for extra security

| Cat. No. | Colour | L x W x H cm | Qty |
|----------|--------|----------------|-----|
| 11780634 | White | 42 x 24 x 20.5 | 1 |

Fisherbrand™ Sharps Container

Safe & Portable

Provides safe disposal of syringes, needles and other small sharps

- 1.0 L capacity with a clear lid to view fill level
- Slide access lid for temporary and final closures for safety and security
- Built-in needle removal port, preventing direct contact with sharps



| Cat. No. | Colour | L x W x H cm | Qty |
|----------|--------|--------------------|-----|
| 15290009 | Yellow | 11.5 x 11.5 x 19.2 | 1 |



Fisherbrand™ Sample Storage Tray

Sample organisation for the busy lab

Hard wearing sample storage trays that can be used for specimen containers with the insert provided, or for general lab products without the insert

- Lightweight and sturdy to hold up to 24 full specimen containers (up to 58 mm in diameter). Specimen containers not included.
- Removable insert to allow outer shell to be used for storage or transportation of laboratory products
- Stackable with other trays to support as many as four high
- Integrated label holder and cards to enable easily changed identification

| Cat. No. | Wells | L x W x H cm | Qty |
|----------|-------|-----------------|-----|
| 15524149 | 24 | 48 x 33.1 x 8.4 | 1 |

Fisherbrand™ Corrugated Sample Storage Tray

Disposable sample organisation and storage

Corrugated storage trays can be used with the included dividers for specimen containers or remove them for general lab product storage or transport

- Made from corrugated polypropylene
- Moisture resistant, significantly reducing the chance of mould

| Cat. No. | Colour | Wells | Well Dia mm | L x W x H cm | Qty |
|----------|--------|-------|-------------|-----------------|-----|
| 15376548 | White | 24 | 58 | 40.6 x 26 x 8.9 | 10 |




Lightweight and Sturdy

to hold up to 24 full specimen containers

Removable Inserts

allow outer shell to be used for storage or transport of laboratory products

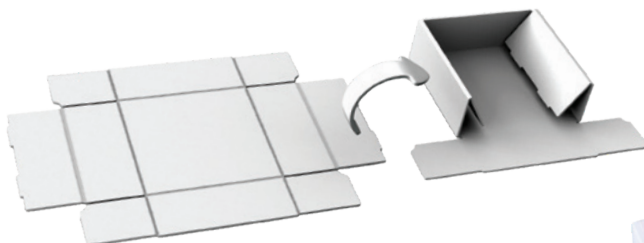

Trays are supplied "flat-packed"

Eliminating the hassle of stacking and reduces the amount of space needed for storage



Carry handle folds
neatly into body space
allowing the units to be
stacked

Moulded Slots
to hold tubes while working



Silicone Water-tight Seal
and three lid clasps provide
secure closure

Fisherbrand™ Transport Box

Specimen or Sample Transport Box

Ideal for transporting biological and clinical samples, as well as instruments and products that you want to keep clean and dry under testing environmental conditions

The transparent box design helps comply with OSHA 29 CFR 1910.1030 standards without having to re-label containers, and enabling the user to visibly check sample integrity.

- Polycarbonate construction and available in three colours for user/owner identification
- Will accommodate 13 and 16 mm tubes in a 72-Place tube rack (page 21) or similar sized/brand of rack
- Use the separators provided to create your own compartment sizes
- Take out separators to transport larger products and equipment
- **Autoclavable** 120°C for 15 min., 16 psi (1.1 bar)



A silicone water-tight seal and three clasps provide secure closure



Load



Lock



Transport



Store




| Cat. No. | Colour | L x W x H cm | Qty |
|----------|---------------------------|-----------------|-----|
| 15165564 | Clear with Blue Handle | ■ | |
| 15175564 | Red with Red Handle | ■ | 1 |
| 15155564 | Yellow with Yellow Handle | ■ | |

Fisherbrand™ Polycarbonate Cryostorage Boxes

Safely store cryogenic vials

Our polycarbonate boxes are designed with unique features to help orient, identify and access vials with ease

- Liquid nitrogen vent and drainage holes
- Magnet on top of forceps enables easy storing on metal surfaces
- Stackable, saves valuable space in the lab
- Store down to temperatures of -196°C with appropriate temperature viable tubes

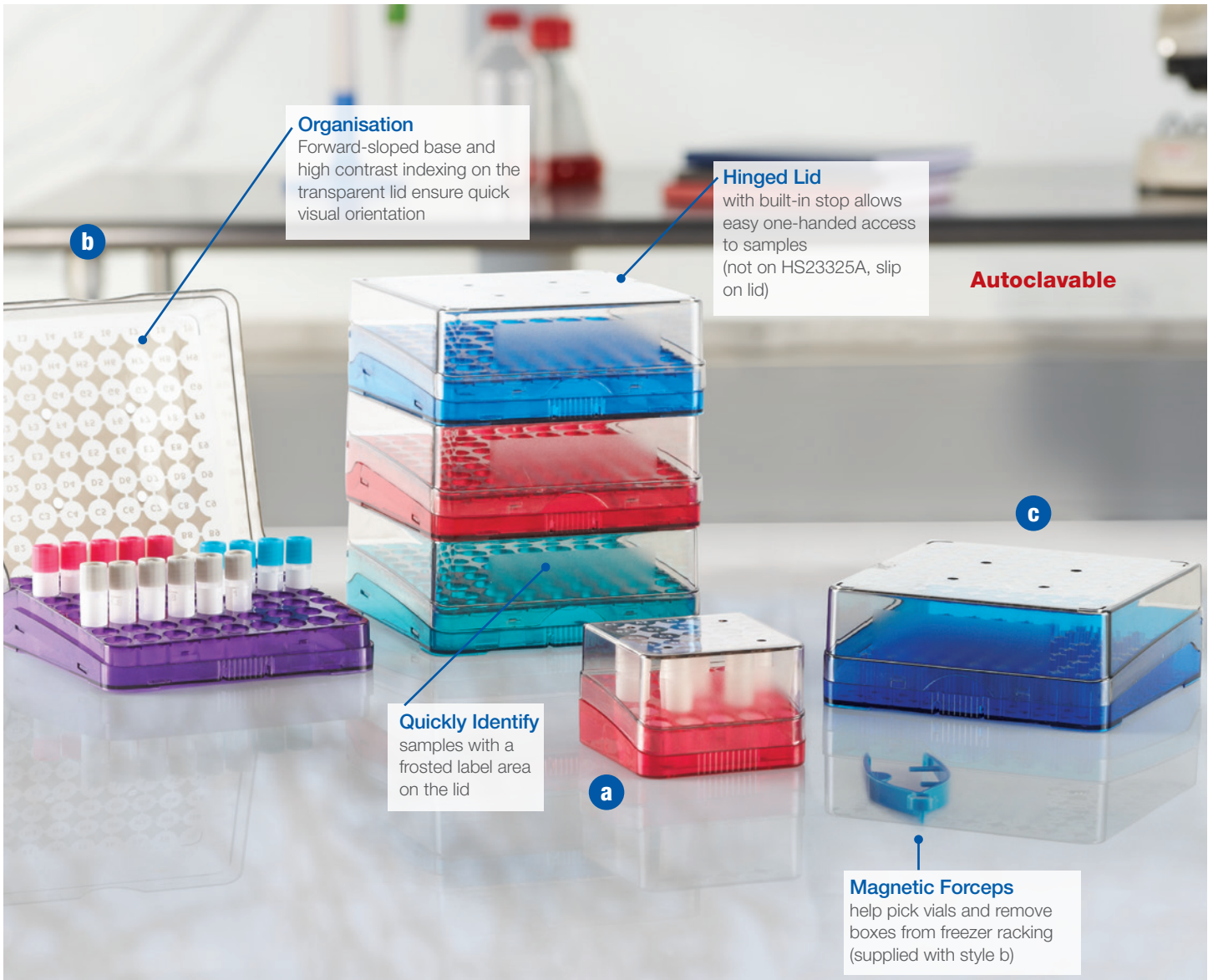

US Patent No. D529622; EP Registered Design No.: 00551429-0001, 00551429-0002, 00551429-0003


| | Cat. No. | Colour | Tube Array | Tubes mL | Wells | Well Dia. mm | Lid | L x W x H cm | Forceps | Qty |
|---|----------|-------------|------------|----------|-------|--------------|----------|-------------------|---------|-----|
| a | 11938084 | Ruby Red | 5 x 5 | 2.0 | 25 | 12.5 | Slide On | 7.6 x 7.6 x 5.3 | | 8 |
| b | 11998004 | Assorted* | 9 x 9 | 2.0 | 81 | 12.5 | Hinged | 13.3 x 13.3 x 5.3 | ✓ | 4 |
| c | 11978004 | Arctic Blue | 10 x 10 | 2.0 | 100 | 12.5 | Hinged | 13.3 x 13.3 x 5.3 | ✓ | 4 |

*Assorted colours (Ruby Red, Amethyst Purple, Arctic Blue, Emerald Green)

Fisherbrand™ Freezer Boxes

Reusable, economical and moisture resistant plastic

- Easy to assemble boxes made from polypropylene material
- Moisture resistant, more durable than cardboard boxes
- Dividers included with the boxes
- Corrugated Polypropylene: Store down to temperatures of -80°C with appropriate temperature viable tubes
- Thin Film: Can be used in mechanical freezers & vessels down to -196°C with appropriate temperature viable tubes. Vapour phase only.

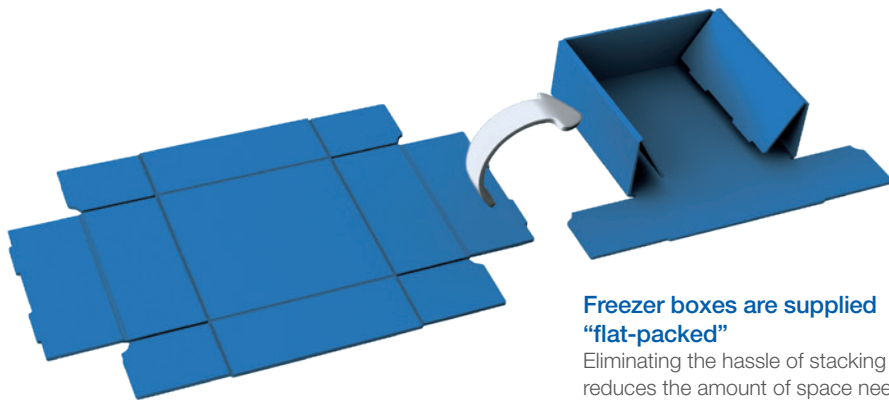

Freezer boxes are supplied "flat-packed"

Eliminating the hassle of stacking and reduces the amount of space needed for storage

Corrugated Polypropylene 0.2/0.5/1.5–2.0/5 mL/Vial



| Cat. No. | Colour | Tubes mL | Wells | Well Dia. mm | L x W x H cm | Qty |
|----------|--------|-----------|-------|--------------|-------------------|-----|
| 15316528 | White | ☐ 0.2 | 144 | 7.9 | 13.3 x 13 x 3 | 10 |
| 15356528 | White | ☐ 0.5 | 81 | 12 | 13.3 x 13 x 4.2 | 10 |
| 15226996 | White | ☐ 1.5/2.0 | 81 | 12 | 13.3 x 13 x 5.1 | 10 |
| 15396528 | Blue | ■ 1.5/2.0 | 81 | 12 | 13.3 x 13 x 5.1 | 10 |
| 15326538 | White | ☐ 5 | 25 | 22.6 | 13.3 x 13 x 7.5 | 10 |
| 15336538 | Blue | ■ 5 | 25 | 22.6 | 13.3 x 13 x 7.5 | 10 |
| 15413790 | White | ☐ Vial | 100 | 17 | 18.8 x 19.8 x 4.1 | 10 |

Corrugated Polypropylene 15 mL and 50 mL

| Cat. No. | Colour | Tubes mL | Wells | Well Dia. mm | L x W x H cm | Qty |
|----------|--------|----------|-------|--------------|------------------|-----|
| 15366538 | White | ☐ 15 | 36 | 19.8 | 14.7 x 14.5 x 12 | 10 |
| 15376538 | Blue | ■ 15 | 36 | 19.8 | 14.7 x 14.5 x 12 | 10 |
| 15306548 | White | ☐ 50 | 16 | 31.6 | 14.7 x 14.5 x 12 | 10 |
| 15316548 | Blue | ■ 50 | 16 | 31.6 | 14.7 x 14.5 x 12 | 10 |

Thin Polypropylene Film 1.5–2 mL

| Cat. No. | Colour | Tubes mL | Wells | Well Dia. mm | L x W x H cm | Qty |
|----------|-----------|-----------|-------|--------------|-------------------|-----|
| 15610785 | White | ☐ 1.5/2.0 | 81 | 12.1 | 12.5 x 12.4 x 4.9 | 10 |
| 15620785 | Blue | ■ 1.5/2.0 | 81 | 12.1 | 12.5 x 12.4 x 4.9 | 10 |
| 15630785 | Lt. Green | ■ 1.5/2.0 | 81 | 12.1 | 12.5 x 12.4 x 4.9 | 10 |
| 15640785 | Purple | ■ 1.5/2.0 | 81 | 12.1 | 12.5 x 12.4 x 4.9 | 10 |

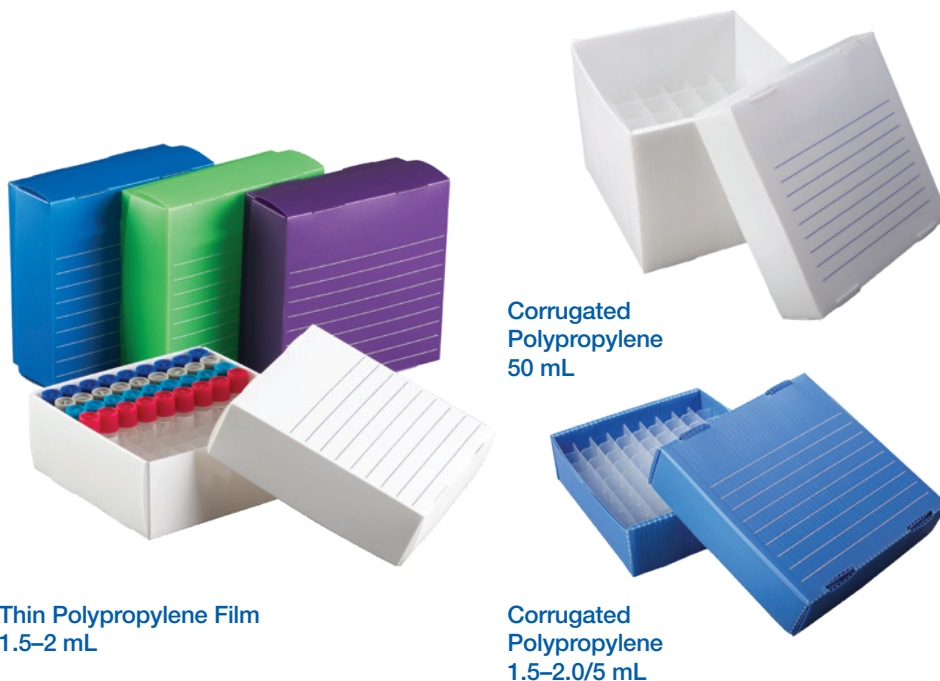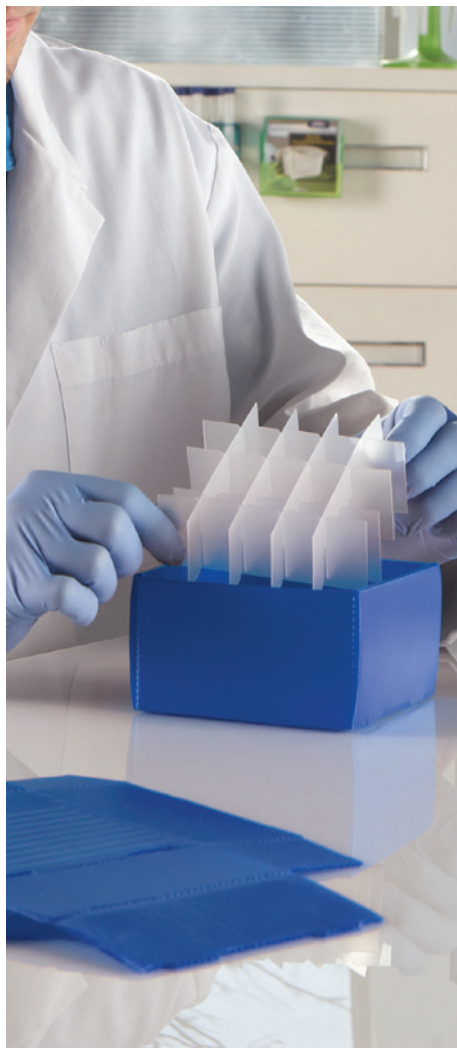

Corrugated Polypropylene 50 mL

Corrugated Polypropylene 1.5–2.0/5 mL

Thin Polypropylene Film 1.5–2 mL

Fisherbrand™ Freezer Microscope Slide Box

Store in ultra-low temperatures with confidence

Protect your slides in low temperature ranges down to -80°C

- Durable polycarbonate and stainless steel clasp/friction fit lid
- Safe storage and transport of 26 x 76 mm slides
- Uniform grooved slide slots separate slides preventing slide contact
- Stay organised with a moulded labelling area with "Project" and "Date"
- Low temperature-resistant numbered foam lining and identification sheet
- Moulded feet enable secure stacking for easy transport and storage



| Cat. No. | Colour | Capacity | Hinge | L x W x H cm | Qty |
|----------|----------|------------|-----------------|-------------------|-----|
| 15366548 | Blue ■ | 25 Slides | Friction Fit | 14.1 x 8.8 x 3.5 | 1 |
| 11730594 | Purple ■ | 100 Slides | Stainless Steel | 20.8 x 17.5 x 3.4 | 1 |


Fisherbrand™ Cryogenic Vial Canes

Organise cryogenic samples in Dewars

Aluminum canes hold vials for freezer storage or immersion in LN₂ Dewars

- Use in liquid and vapour phase; use approved vials only
- Canes bend to fit all major cryogenic vial brands

| Cat. No. | Tubes mL | No. of Vials | Vial Size | Length cm | Qty |
|----------|----------|--------------|-----------|-----------|-----|
| 11958084 | 15 & 50 | 5 | Up to 2.0 | 2.9 | 12 |
| 11968084 | 15 & 50 | 6 | Up to 5.0 | 3 | 12 |

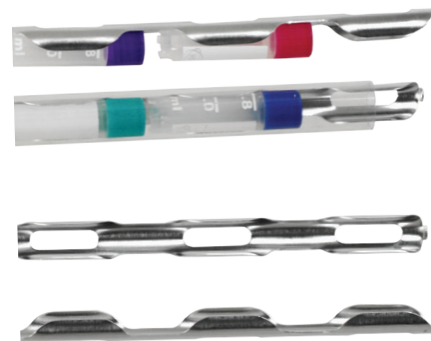

Fisherbrand™ PVC Cryogenic Cane Sleeves

For extra security and sample safety

Clear plastic sleeves protect vials in aluminum canes during handling and storage

- Safe for cryogenic freezing
- More durable than cardboard sleeves
- Transparency allows sample identification without removing sleeve

| Cat. No. | Diameter cm | Length cm | Qty |
|----------|-------------|-----------|-----|
| 11734673 | 1.5 | 27.5 | 100 |



Fisherbrand™ Cryogenic & Vaccine Vial Grippers

These two styles of hand-operated grippers are designed to pick up cryogenic vials and crimped vaccine/serum vials by the cap, respectively.

Made to grab with precision and safely remove tubes and vials from cold temperature storage boxes.

Both grippers are crafted from autoclavable, reinforced polypropylene, and feature an ergonomic squeeze area for optimal precision.

Each style is available in a pack of 5 assorted colors for easy identification and organization. A combo pack with one of each style is also available.

Autoclavable 120°C for 15 min



| Cat. No. | Type | Colour | L x W x H cm | Qty |
|----------|-----------|---------------------|------------------|-----|
| 17807155 | Cryogenic | Assorted* ■ ■ ■ ■ ■ | 14.3 x 1.5 x 2.3 | 5 |
| 17887145 | Vaccine | Assorted* ■ ■ ■ ■ ■ | 15.2 x 2 x 3.5 | 5 |

Assorted colors: red, blue, green, purple, yellow

| Cat. No. | Type | Colour | L x W x H cm | Qty |
|----------|-------------|--------|--------------------|-----|
| 17897145 | Combination | Blue ■ | 45.9 x 15.9 x 27.1 | 2 |

Fisherbrand™ Universal Freezer Rack

Innovative horizontal drawer system for upright freezers

- Constructed of corrosion-resistant stainless steel with four drawers to hold 16 boxes
- Drawers feature a unique stop, that when fully extended pivot down allowing for easy box access
- Steel drawer handles provide a user-friendly grip
- Freezer rack also includes a pull-out wire handle for easy removal off the freezer shelf
- Compatible with most standard 5.3 cm freezer storage boxes



Unique drawer stop that pivots down allowing for easy box access


Steel drawer handles have a space to insert a label, and provide a user-friendly grip

| Cat. No. | Config. | Box Height cm | L x W x H cm | Qty |
|----------|---------|---------------|--------------------|-----|
| 11786013 | 4 x 4 | 5.3 | 13.9 x 57.1 x 23.8 | 1 |

Fisherbrand™ Cube Freezer Rack – Upright

For 15 and 50 mL tube box storage

- Constructed of corrosion-resistant stainless steel with 6 rack spaces
- Provides safe, convenient storage of your samples
- Handles at each end accommodate freezers with right or left handed doors
- Compatible with 15.2 cm or smaller boxes

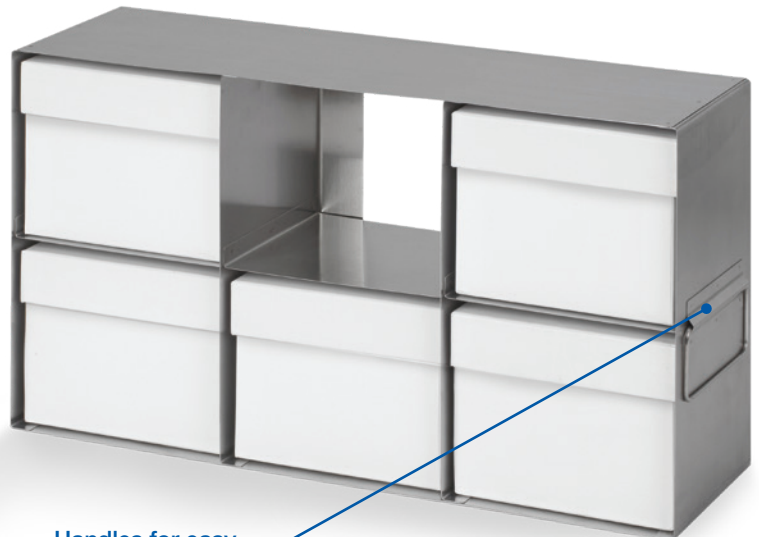

Handles for easy transfer

Robust mounting brackets assure handles stay securely in place for easy use

| Cat. No. | Config. | Opening L x W x H cm | L x W x H cm | Qty |
|----------|---------|----------------------|--------------------|-----|
| 12643476 | 2 x 3 | 15.2 x 15.7 x 13.4 | 45.9 x 15.9 x 27.1 | 1 |



Fisherbrand™ Paper Boats

Green solution to sample weighing

Made from FSC™ Certified (FSC™ Mix Credit) paper (Eco-IN) and when disposed of in the standard waste stream they will naturally biodegrade (Eco-OUT)

- A cost-effective solution to reduce plastics in the lab
- Made of grease-proof paper that meets FDA 21 CFR 176.170 (components of paper in contact with aqueous and fatty foods) and with 21 CFR 176.180 (contact with dry foods)
- Control contents like a polystyrene weigh boat because they are easier to pour
- An alternative to regular weighing papers because they are easier to handle

| Cat. No. | Colour | Size | L x W x H mm | Qty |
|----------|--------|--------|--------------|-----|
| 15422214 | Brown | Small | 35 x 35 x 15 | 500 |
| 15408113 | White | Small | 35 x 35 x 15 | 500 |
| 15433214 | Brown | Medium | 45 x 45 x 15 | 500 |
| 15442214 | White | Medium | 45 x 45 x 15 | 500 |
| 15893967 | Brown | Large | 95 x 95 x 25 | 250 |
| 15853967 | White | Large | 95 x 95 x 25 | 250 |

Functional in the Lab

Available in a small and medium in both brown for light coloured samples and white for dark samples



Fisherbrand™ Disposable Spatulas


Designed to fit perfectly with standard 0.2, 0.5, 1.5/2.0, 15 and 50 mL tubes

- Made of virgin, food-grade polypropylene
- Ideal for dispensing or sample collection of powders or soft solids
- Durable to work with liquid nitrogen frozen samples
- Hollow core minimises heat transfer; good for hot and cold liquids
- Anti-static, sterile and eco-friendly versions available

Always there and always clean

Autoclavable

Lab tools that scoop, stir and scrape

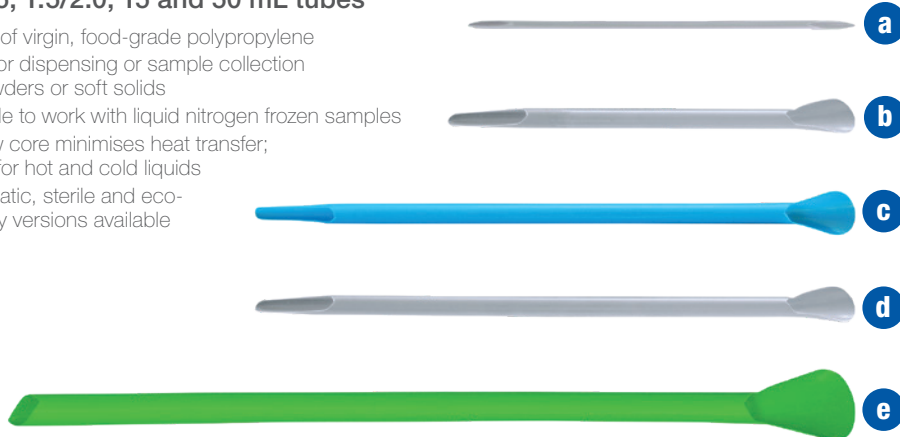


| Cat. No. | Description | Colour | Suggested Tube mL | L x W cm | Spoon (Flare) mm | Scoop/Pick L x W mm | Taper | Qty |
|------------|-------------------|--------|-------------------|----------|------------------|---------------------|-------|-----|
| a 15356558 | Micro Anti-static | Opaque | 0.2 and 0.5 | 14 x .4 | — | 11 (V-Tip) | 15° | 300 |
| b 15326134 | Eco-Friendly | Opaque | 15 and 1.5 | 14 x .7 | 25 x 14 | 30 x 3 | | 300 |
| c 15336558 | Standard | Blue | | | | | | 300 |
| d 15346558 | Standard | Opaque | 15 and 1.5 | 21 x .7 | 25 x 14 | 30 x 3 | | 300 |
| 15376558 | Std. EtO Sterile | Opaque | | | | | | 100 |
| e 15366558 | Macro | Green | 50 | 31 x 1 | 42 x 23.5 | 30 | 35° | 150 |

Fisherbrand™ Weighing Boats

Flexible polystyrene boats feature smooth surfaces to provide accurate pour-outs with minimal sample loss

- Ideal for weighing, dispensing or storing applications
- Select white, black or blue for best colour contrast to ensure sample is completely removed
- Biologically inert and resistant to dilute and weak acids, aqueous solutions, alcohols and bases
- Temperature range -10 to 70°C (14 to 158°F)
- Antistatic versions available making the boats non-reactive with most substances and suitable for weighing static-affected samples
- Wide flat bottoms ensure stability and prevent tipping
- Rounded corners allow simple flexing for controlled pouring of liquid, powdered or granular samples
- Flat bottoms and sloping sides allow easy stacking for convenient storage
- Can be used with microwave moisture analysers
- Made from FDA/food-grade compliant material

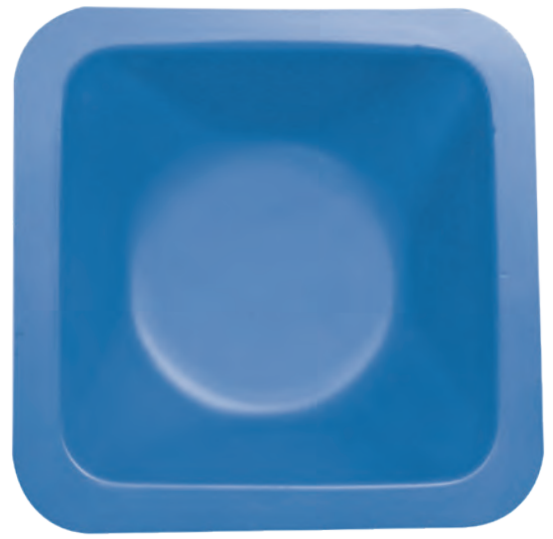

Standard Weighing Boats



Sterile Weighing Boats

Ideal for applications where there is concern of cross contamination and 'clean' is a necessity in life science research, food, bio-pharmaceutical, environmental and molecular research environments

- Electron beam sterilized
- Lot number traceability for quality assurance
- Certificate of Sterility available upon request

STERILE R

NEW



Standard




Hexagonal



Pour-boat



Square

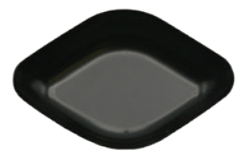

Diamond

Weigh Boats

| Shape | Regular Material Cat. No. | Antistatic Material Cat. No. | Colour | Size | Volume mL | L x W x H cm | Qty |
|-----------|---------------------------|-------------------------------------------|-------------------------------------------|---------|-----------|-------------------|-----|
| Standard | 11593422 | 12922850 | White <input type="checkbox"/> | Small | 10 | 4.6 x 4.6 x .8 | 500 |
| | 15728187 | 15768167 | Blue <input checked="" type="checkbox"/> | | | | |
| | 12982850 | 12942860 | Black <input checked="" type="checkbox"/> | | | | |
| | 11503432 | 11680302 | White <input type="checkbox"/> | | | | |
| | 15738187 | 15778167 | Blue <input checked="" type="checkbox"/> | Medium | 85 | 8.5 x 8.5 x 2.4 | 500 |
| | 12992850 | 12952860 | Black <input checked="" type="checkbox"/> | | | | |
| | 11513432 | 12608513 | White <input type="checkbox"/> | Large | 280 | 14 x 14 x 2.2 | 500 |
| | 15748187 | 15788167 | Blue <input checked="" type="checkbox"/> | | | | |
| 12902860 | 12962860 | Black <input checked="" type="checkbox"/> | | | | | |
| Hexagonal | — | 15778187 | White <input type="checkbox"/> | Small | 6 | 3.8 x 2.5 x .9 | 500 |
| | — | 15788187 | | Medium | 46 | 6.4 x 5.1 x 1.9 | 500 |
| | — | 15798187 | | Large | 194 | 10.2 x 7.6 x 2.2 | 500 |
| | — | 15708197 | | X-Large | 364 | 12.7 x 8.9 x 3.1 | 500 |
| Pour-boat | 15798177 | 15798167 | White <input type="checkbox"/> | Small | 14 | 4.3 x 5.8 x 1.3 | 500 |
| | 15708187 | 15708177 | | Medium | 137 | 8.3 x 13.2 x 2.6 | 500 |
| | 15718187 | 15718177 | | Large | 259 | 10.8 x 18.3 x 2.6 | 500 |
| Square | 15758187 | — | White <input type="checkbox"/> | Medium | 51 | 9 x 9 x .7 | 500 |
| Diamond | 12932840 | 12992840 | White <input type="checkbox"/> | Small | 5 | 5.3 x 3.3 x .6 | 500 |
| | 12952850 | 12912860 | Black <input checked="" type="checkbox"/> | | | | |
| | 15768187 | 15718197 | White <input type="checkbox"/> | Medium | 25 | 6.7 x 4.7 x 1.4 | 500 |
| | 15738197 | 15728197 | Black <input checked="" type="checkbox"/> | | | | |
| | 11573422 | 11533432 | White <input type="checkbox"/> | Large | 30 | 7.8 x 5.6 x 1.4 | 500 |
| | 12962850 | 12922860 | Black <input checked="" type="checkbox"/> | | | | |
| | 12387552 | 11543432 | White <input type="checkbox"/> | X-Large | 100 | 11.9 x 9 x 1.9 | 500 |
| | 12972850 | 12932860 | Black <input checked="" type="checkbox"/> | | | | |

Sterile Weigh Boats


| Shape | Antistatic Material Cat. No. | Color | Size | Volume mL | L x W x H in | L x W x H cm | Qty |
|-----------|------------------------------|--------------------------------|--------|-----------|----------------|-----------------|----------------------|
| Standard | 17754824 | White <input type="checkbox"/> | Small | 10 | 1.8 x 1.8 x .3 | 4.6 x 4.6 x .8 | Pack of 200/40 x 5pk |
| Hexagonal | 17764824 | White <input type="checkbox"/> | Small | 6 | 1.5 x 1 x .4 | 3.8 x 2.5 x .9 | |
| Hexagonal | 17774824 | White <input type="checkbox"/> | Medium | 46 | 2.5 x 2 x .7 | 6.4 x 5.1 x 1.9 | |

Fisherbrand™ Liquid Pour Boats

Deep vessel with sturdy handle designed for convenient handling and pouring of liquids

- Flat bottom ensures stability and prevents tipping
- Notched pour spout design controls liquid pour
- Operating temperature range 0 to 40°C (104°F)
- Made from FDA food grade polystyrene in an ISO Class 8 clean room

- Material biologically inert and resistant to dilute and weak acids, aqueous solutions, alcohols and bases
- Recyclable, marked
- Available in 3 Sizes
 - 50mL
 - 100mL
 - 200 mL

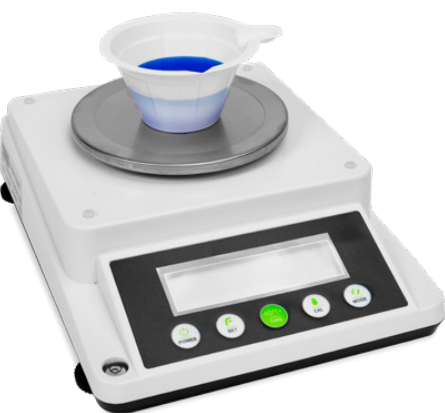

NEW



Built-in Pour Spout
for accurate pouring

Large Sturdy Handle
for ease of movement and pouring

| Cat. No. | Color | Size | Volume mL | L x W x H in | L x W x H cm | Qty |
|----------|--------------------------------|--------|-----------|-----------------|-------------------|-------------|
| 17860752 | White <input type="checkbox"/> | Small | 50 | 3.6 x 1.6 x 1.7 | 9.2 x 4.0 x 4.3 | Pack of 200 |
| 17870752 | White <input type="checkbox"/> | Medium | 100 | 4.2 x 3.5 x 1.9 | 10.7 x 9 x 4.8 | |
| 17880752 | White <input type="checkbox"/> | Large | 200 | 5.8 x 4.6 x 2 | 14.8 x 11.7 x 5.1 | |

Fisherbrand™ Micro Weigh Boats

NEW

Ideal for microbalance weighing and/or distribution of very small amounts of sampling material


Easy to Handle

- Large tabs make boats easy to use, position and hold, especially with forceps

Superb Sample Recovery

- Inert material enables the boat and contents to be added to solutions to aid 100% sample recovery, where protocols allow
- Anti-static properties prevent sample retention from static electricity
- Small surface area reduces liquid sample evaporation

| Item No. | Color | Dimensions mm | Qty |
|----------|-------|----------------|-----|
| 17113855 | White | 25 x 8.9 x 7.5 | 250 |

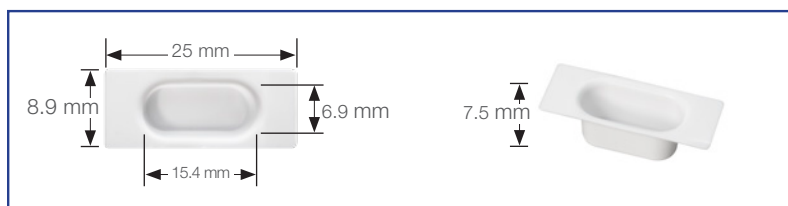

Actual size

Easy to handle tabs



Holds approx 500 µl

Made of polystyrene
Food Grade Compliant



Ideal for very small volume pipette calibration and checks and micro sampling



NEW

Fisherbrand™ Dual Purpose Scoops

Scoop + Weigh Boat

Unique economical solution that combines the function of a scoop and a weigh boat into one, streamlining the process. Just scoop to retrieve, place on scale to weigh

- Made from anti-static polystyrene that is waterproof and reduces the chance of materials sticking to the scoop
- Ships and stores flat, easy to assemble for use
- Locking tab feature enables a scoop shape to help retain materials
- Coated to stand up to powders, gels and semi-solids
- Electronic Beam Sterilisation procedures compliant with ISO11137-3:2017
- Certificate of Sterility available upon request



| Cat. No. | Colour | Size | L x W cm | Qty | |
|----------|--------|--------------------------|----------|------------|-----|
| 15833977 | White | <input type="checkbox"/> | Large | 12.7 x 9 | 100 |
| 15843977 | White | <input type="checkbox"/> | Small | 12.7 x 6.3 | 100 |

Fisherbrand™ Weighing Boats & Dishes

For general-purpose weighing, dispensing, storage, food processing, evaporating or drying

Aluminum Weighing Boats

- Sturdy yet lightweight boats feature a smooth, flat bottom for sample handling and a flanged, crimped wall for rigidity
- Finger tab handle can be etched for identification prior to weighing



Aluminum Weighing Boats




Aluminum Weighing Dish

Aluminum Weighing Dish

- Designed for moisture determination balances with slightly tapered walls for stacking
- Use as dust covers for open vessels



| Cat. No. | Size | Dia. x H mm | Qty |
|----------|--------|-------------|-----|
| 11998084 | Small | 43 x 12 | |
| 11908094 | Medium | 57 x 14 | 100 |
| 11918094 | Large | 70 x 16 | |
| 11988084 | Large | 10 x 8 | 50 |

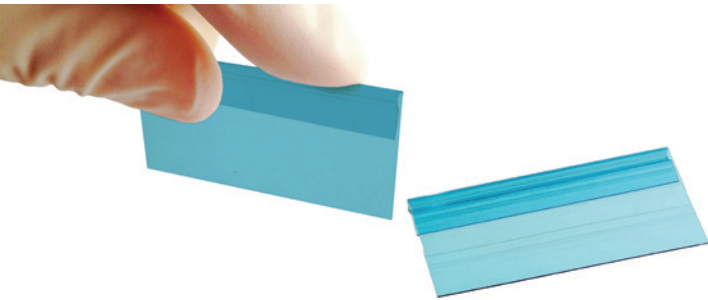


Fisherbrand™ Safe Slicer

Laboratory Safety Razor Blade

Say no to dangerous sharps

- At just 0.27 mm the US Patent Pending product is as thin as a single edge industrial razor blade, but much safer
- Cleanly slice directly on light box without scratching the surface
- Adjust the size required by simply cutting with scissors
- Made from food grade PET

| Cat. No. | L x W x H mm | Qty |
|----------|---------------|-----|
| 15316558 | 19 x 38 x 0.3 | 100 |
| 15326558 | 25 x 50 x 0.3 | |

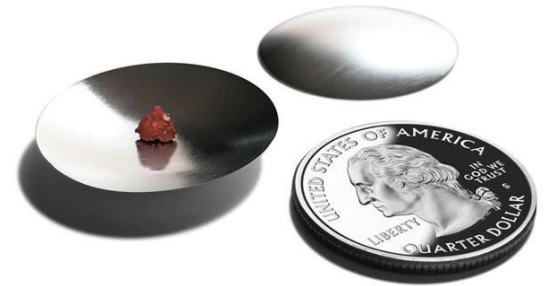

Fisherbrand™ Aluminum Weighing Dish

- Ideal for microbalances, organic micro elementary and thermal analysis, and pharmaceutical use
- Made of Ultra High Vacuum (UHV) Foil - "A" wettable and completely free from oils normally found on other grades of aluminum foil
- Foil meets ASTM B-479, specification 3.1.4 and 10.3.1
- "Cahn style" accessory pans are compatible with most commercially available analysers and microbalances

| Cat. No. | Dia. x H mm | Qty |
|----------|-------------|-----|
| 15803977 | 25 x 3.6 | 20 |

Fisherbrand™ Paper Disposable Funnel

A clean transfer every time

- Supplied flat-pack to minimise storage
- With a minimal amount of finger pressure the funnel takes its working shape
- Maximum flexibility – cut the stem to any size or shape to meet with your needs
- Fits a wide variety of openings and can be used with virtually any liquid* or powder
- *Not recommended for hot liquids
- Disposable



| Cat. No. | Capacity mL | Top Dia. cm ² | Bottom Dia. cm ² | Height cm | Qty |
|----------|-------------|--------------------------|-----------------------------|-----------|-----|
| 15346568 | 186 | 6 | 0.6 | 12.8 | 25 |

Fisherbrand™ Lead Rings

Stabilise vessels and labware or immerse in water baths

- Vinyl-coated lead rings will not mark or scratch benchtops and labware
- C-shape format allows some manual size variation for a customised fit



| Cat. No. | Shape | Colour | Flask Sizes mL | ID cm | OD cm | Weight g | Qty |
|----------|-------|--------|----------------|-------|-------|----------|-----|
| 11726725 | C | Blue | 125–500 | 4.2 | 7.5 | 215 | 1 |
| 10195573 | O | Blue | | 4.8 | 7.7 | 210 | 1 |
| 11736725 | C | Blue | 250–1000 | 5.1 | 9 | 465 | 1 |
| 11797884 | O | Blue | | 5.1 | 9 | 490 | 1 |
| 11746725 | C | Blue | 500–2000 | 5.4 | 9.8 | 595 | 1 |
| 11707894 | O | Blue | | 5.7 | 9.8 | 610 | 1 |
| 11756725 | C | Blue | 1000–4000 | 6.6 | 11.3 | 890 | 1 |
| 11717894 | O | Blue | | 7 | 11.5 | 935 | 1 |

*Weights may vary by +/- 3%

Fisherbrand™ Cell Spreader

- Made from ABS
- Flexible rippled handle and domed spreading bar to help reduce agar tears
- Gamma sterilised

Domed Headpiece
promotes smooth spreading




| Cat. No. | Shape | L x W x H cm | Packaging | Sterile | Qty |
|----------|-------|-----------------|------------|---------|-----|
| 11364035 | Delta | 17.2 x 3.7 x .4 | 1 each/bag | ✓ | 500 |

Fisherbrand™ Inoculating Picks

The 'smart' way to isolate individual colonies


Available in three lengths and two needle sizes to suit a wide range of application needs

- Naturally biodegradable, made from 100% U.S.A. grown white birch wood
- Manufactured in an FDA-registered factory and complies with FDA standards for food contact
- Available in semi or needle point styles for delicate operations

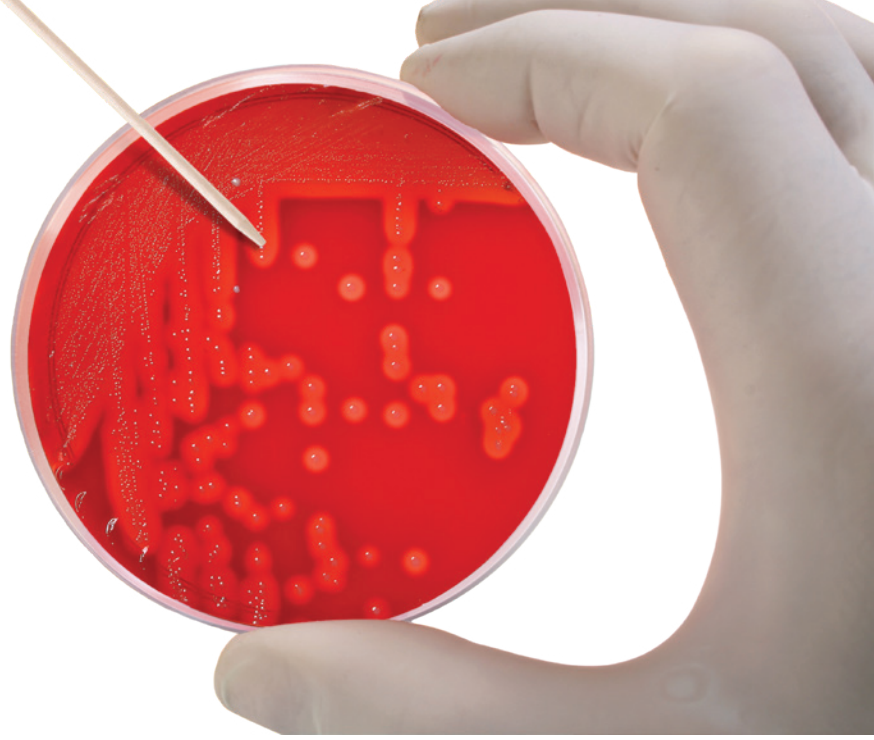

Autoclavable

| Cat. No. | Size | Point | Tip | Suggested Tube mL | Length cm | Qty |
|----------|----------|--------------|--------|----------------------|-----------|------|
| 15396558 | Mini | Needle Point | | Fits 1.5 Closed Tube | 3.3 | 2000 |
| 15306568 | Standard | Semi Point | Single | | 6.4 | 1600 |
| 15316568 | Standard | Needle Point | Tip | | | |

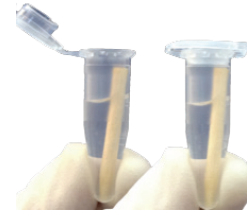

Mini Sized Inoculating Picks

fit neatly into a 1.5 mL microfuge tube allowing the lid to be closed for full sample transfer/recovery

Drop it in & Close it up



Fisherbrand™ Workstation Dispensing Bins

Crystal-clear bins let you see your small items inventory

- See inventory levels of pipette tips, tubes, fittings, finger cots and other bulk small items at a glance
- Flip-top covers on generous openings make filling and dispensing easy and keep products dust-free
- 4 sizes available to cater to a wide variety of storage needs
- Includes non-skid rubber feet

| Cat. No. | Colour | Capacity | L x W x H cm | Qty |
|----------|--------------------------------|----------|--------------------|-----|
| 11764834 | Clear <input type="checkbox"/> | 1 Bin | 18 x 17.5 x 23.5 | 1 |
| 11774834 | | 2 Bins | 13.9 x 16 x 36.8 | |
| 11907994 | | 3 Bins | 21 x 15.5 x 21.2 | |
| 11784834 | | 4 Bins | 27.8 x 15.8 x 21.7 | |

Fisherbrand™ Petri Dish Dispenser

- Sturdy acrylic holds up to 30 x 100 mm dishes in two easy-access partitions
- Bench top or wall mountable (screws provided)



| Cat. No. | Colour | L x W x H cm | Qty |
|----------|--------------------------------|--------------------|-----|
| 11384145 | Clear <input type="checkbox"/> | 22.2 x 10.7 x 23.6 | 1 |

Fisherbrand™ Lab Wipe Holders

Holds Lab Wipes such as Kimwipes™ and other major brands
Clear acrylic holders fit boxes of small or large lab wipes.

- Mount on a shelf or wall
- The smaller version attaches with double-stick tape and the larger version with screws
- Mounting materials included



| Cat. No. | Size | L x W x H cm | Qty |
|----------|-------|------------------|-----|
| 12377619 | Small | 1.8 x 12.2 x 8.4 | 1 |
| 12387619 | Large | 32 x 12.7 x 9.7 | |

Fisherbrand™ Clear Acrylic Glove Box Holders

Keep gloves organised & ready for immediate use
Free up bench-top space by mounting holders on wall for easy dispensing of gloves

- Universal fit for most major brands of glove boxes
- Mounting screws included
- Includes a polyester foam insert to securely hold a range of glove box sizes—snug fit ensures easy removal of gloves from the box

| Cat. No. | Capacity | L x W x H cm | Qty |
|----------|----------|--------------------|-----|
| 11897102 | 1 box | 26.3 x 14.6 x 11.4 | 1 |
| 11807112 | 2 boxes | 29.2 x 26 x 11.4 | 1 |
| 11817112 | 3 boxes | 29.4 x 38.8 x 11.6 | 1 |



Fisherbrand™ Epoxy Coated Glove Box Holders

Keep gloves organised & ready for immediate use

- Epoxy-coated steel wire racks are an economical solution
- Universal fit for most major brands of glove boxes
- Mounting screws included

| Cat. No. | Capacity | L x W x H cm | Qty |
|----------|----------|------------------|-----|
| 11885843 | 1 box | 14 x 20.9 x 10.8 | 1 |
| 11895843 | 3 boxes | 25 x 45.4 x 10.8 | 1 |

Fisherbrand™ Modular Glove Box Holders

Ideal for use in clean rooms, hospitals, and microbiological work areas.

- High-impact, durable resin is resistant to harsh cleaning solutions
- Modular design allows single boxes or horizontal stacking for multiple box holders
- Universal fit for most major brands of glove boxes
- Mounting screws included

| Cat. No. | Description | Colour | L x W x H cm | Qty |
|----------|----------------|------------------------------------------|-------------------|-----|
| 15380440 | Standard | White <input type="checkbox"/> | 25.4 x 15.6 x 9.7 | 3 |
| 11710644 | Anti-Microbial | Blue <input checked="" type="checkbox"/> | | 3 |



US Patent No. 6669157



Adaptable

Store different size boxes all in one holder

Fisherbrand™ Flexible Glove Box Holder

Glove box holder ideal for gloves, face masks, shoe covers...

The adjustable holding points provide true flexibility to hold up to four different boxed protective products and other lab essentials (size dependent)

- Ideal for holding different boxed sizes of gloves, mop caps, ear plugs, face masks, shoe covers and more
- Adjustable universal fit
- Easy-to-clean, with a low profile
- Wall- or door-mountable
- Mounting hardware included

| Cat. No. | Colour | L x W x D cm | Qty |
|----------|--------------------------------|-------------------|-----|
| 15279568 | White <input type="checkbox"/> | 56.5 x 13.3 x 4.6 | 1 |



Fisherbrand™ Workstation Storage Bin

The essential combination accessory holder

Organise your most commonly used laboratory disposables in one central location

- Acrylic unit has room to store gloves, Kimwipes™ and Parafilm™
- Use bins to store tubes, tips, etc.

| Cat. No. | Colour | L x W x H cm | Qty |
|----------|--------------------------------|------------------|-----|
| 11927994 | Clear <input type="checkbox"/> | 26.5 x 13.7 x 44 | 1 |

Parafilm™ is a registered trademark of Bemis Inc.

Fisherbrand™ Drawer Organiser

Keeps supplies organised and easily accessible

Polystyrene organiser perfect for storing and organising small laboratory supplies

- Protect small and fragile items
- Fit in most standard drawers or store on benches and shelves



| Cat. No. | Colour | L x W x H cm | Qty |
|----------|--------|-------------------|-----|
| 11747444 | White | 45.5 x 30.7 x 6.3 | 1 |

| Compartments | L x W cm |
|----------------------|----------|
| 2 Small Compartments | 43 x 8 |
| 1 Large Compartment | 43 x 11 |

Fisherbrand™ Storage Racks for 80-Well Racks

Strong, clear acrylic construction

These racks help you keep large numbers of microtubes organised on your benchtop

- Holds up to eight 80-well microtube racks filled with either 1.5 or 2.0 mL microtubes
- Plenty of room to slide the racks in and out



| Cat. No. | Description | Colour | L x W x H cm | Qty |
|----------|-------------|--------|--------------------|-----|
| 11708214 | Vertical | Clear | 15.7 x 23.6 x 23.4 | 1 |
| 11710055 | Horizontal | Clear | 31 x 23.6 x 12.2 | |

Fisherbrand™ Benchtop Dispensing Bins

Work great with gloved hands, bright colours

Ergonomically designed with large, angled openings and oversized lip

- Lid remains open when tilted back
- Non-skid rubber feet and pre-drilled holes for wall mounting
- Green bin includes a removable divider
- Mounting screws included



| Cat. No. | Colour | Compartments | L x W x H cm | Qty |
|----------|------------|--------------|------------------|-----|
| 11304065 | Neon Green | 2 | 18 x 15.5 x 16.5 | 1 |
| 11314065 | Neon Blue | 1 | 15.5 x 15.5 x 17 | |

Fisherbrand™ Earplug Dispenser

Translucent acrylic dispensers provide easy access to safety accessories

- Holds approximately 200 pairs of foam earplugs
- Set on counter or wall mount (screws included)

| Cat. No. | Colour | L x W x D cm | Qty |
|----------|--------------------------------------------|--------------------|-----|
| 11507473 | Green ■ | 22.7 x 16.2 x 15.6 | 1 |



Fisherbrand™ Soft Cover Dispenser

Store and dispense with ease

Translucent acrylic dispensers provide easy access to safety accessories

- Hinged lid protects contents from dust and debris
- Set on counter or wall mount (screws included)
- Great for holding other personal items

| Cat. No. | Colour | L x W x D cm | Qty |
|----------|---------------------------------------------|--------------------|-----|
| 15599801 | Natural □ | 22.7 x 16.2 x 15.6 | 1 |

Fisherbrand™ 3-in-1 Dispenser


Three separate lift-up lids make loading and dispensing supplies easy

Three-compartment dispenser can be used to organise safety glasses, earplugs, soft covers (including bouffant caps, beard covers, shoe covers) or any small items that need to be kept within reach

The unit includes screws for wall mounting or can be set on a counter

- Middle compartment has reach-in hole for dispensing
- Durable acrylic is see-through so product levels are always visible

| Cat. No. | Colour | L x W x D cm | Qty |
|----------|------------------------------------------|--------------------|-----|
| 11927984 | Gray ■ | 40.6 x 40.6 x 20.3 | 1 |

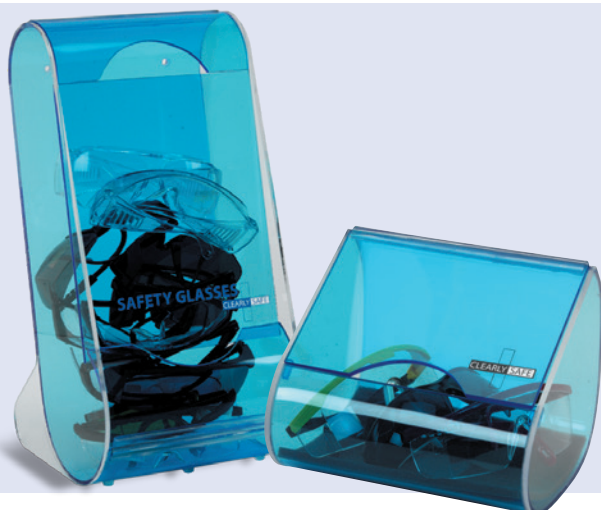

Fisherbrand™ Safety Glasses Dispensers

Ergonomically angled for easy dispensing

Translucent acrylic dispensers provide easy access to safety accessories

- Hinged lid protects glasses from dust and debris
- Set on counter or wall mount (screws included)
- Large dispenser holds approx. 20 glasses, small approx. 10 glasses
- Great for holding other personal items

| Cat. No. | Colour | L x W x D cm | Qty |
|----------|------------------------------------------|--------------------|-----|
| 11803470 | Blue ■ | 22.7 x 16.2 x 15.6 | 1 |
| 11893460 | Blue ■ | 20.3 x 20.3 x 40.6 | 1 |



Constructed of high-impact polypropylene

Fisherbrand™ Handy Container

Sturdy storage tub with lid and divider

Our handy stackable storage container is designed to help organise a laboratory

- The removable dividers can create as big or as small a compartment as needed. Perfect for keeping tubes, tips, pipettes, gloves and other small items organised and out of the way
- The lid snaps in place to secure the contents and for safe transport

| Cat. No. | L x W x H cm | Qty |
|----------|----------------|-----|
| 11938014 | 33 x 20 x 11.5 | 5 |



Fisherbrand™ MagLab

Magnetic lab storage series

Easily organise common lab products off your bench top with the MagLab storage system

Vibrant & Functional

- Made of robust ABS plastic in bright translucent colours
- Secure direct to metal with built-in powerful magnets or non-metal surfaces with included self-adhesive mounting plate
- Available individually or as a pack of four

| Cat. No. | Description | Colour | L x W x H cm | Qty |
|----------|---------------|----------|-------------------|-----|
| 11942581 | MagLab Series | Assorted | — | 4 |
| 11992581 | MagWipe | Green | 13 x 9.1 x 9.8 | — |
| 11932581 | MagPette | Purple | 13.1 x 11.2 x 7.6 | 1 |
| 11922591 | MagBox | Orange | 13 x 7.8 x 10.2 | — |


Kimwipes™ is a registered trademark of Kimberly-Clark Corp. *Assorted colours: Teal, Green, Purple and Orange

MagRack holds four 50 mL tubes and nine 15 mL tubes (only available in MagLab Series)

MagWipe holds a standard (110 x 210 mm) Kimwipes™ box

MagPette holds two pipettors (fits most popular styles)

MagBox holds miscellaneous accessories

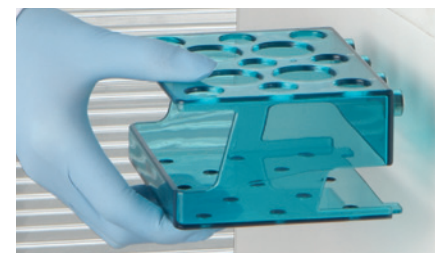

MagRack (only available in MagLab Series)

Fisherbrand™ Box Top Dispenser

Easy application with no mess

Easy application, simply place the dispenser on top of the Parafilm® M box and pull the sealing film through, there is no need to remove rolls from the packaging

- Accommodates two 50 mm or one 100 mm Parafilm® M rolls
- Built in safety blade slides and cuts with no danger to fingertips and eliminates the need for scissors
- Supplied with two safety blades, additional replacement blades are available for purchase
- Made from ABS plastic making the unit easy to clean

NEW



| Cat. No. | | L x W x H cm | Qty |
|----------|--------------------------------------------------------------------------------------|-------------------|-----|
| 16390412 | Natural <input type="checkbox"/> | 12.3 x 17.2 x 9.9 | 1 |
| 16390402 | Replacement blades for Box Top Dispenser - Black <input checked="" type="checkbox"/> | | 3 |

Parafilm® is a registered trademark of Amcor PLC.



Fisherbrand™ ABS dispenser for Parafilm™ Sealing Film

Unique design for maximum safety

Store, dispense and cut Parafilm™, tape, labels and Tough-Tags™

- ABS plastic dispenser designed with a safety razor that slides and cuts with no danger to fingers and supplied with two easily exchanged cutting blades
- Dispenser holds one or two 50 mm rolls or one 100 mm roll

| Cat. No. | | L x W x H cm | Qty |
|----------|------------------------------------------|------------------|-----|
| 11350040 | Blue <input checked="" type="checkbox"/> | 12 x 15.6 x 17.1 | 1 |
| 11865993 | Green <input type="checkbox"/> | | |
| 11875993 | Red <input type="checkbox"/> | | |
| 11928014 | Replacement blades for ABS dispenser | | 5 |

Parafilm™ is a registered trademark of Amcor PLC. Tough-Tags™ is a registered trademark of Diversified Biotech, Inc.

Fisherbrand™ Acrylic Dispenser for Parafilm™ Sealing Film

Heavy-duty construction for increased durability

Store, dispense and cut Parafilm™, tape, label, and Tough-Tags™

- Strong, crystal-clear acrylic dispenser makes your film supply visible and features a serrated cutting blade
- Dispenser holds one or two 50 mm (2-inch) rolls or one 100 mm (4-inch) roll

| Cat. No. | L x W x H cm | Qty |
|----------|--------------------|-----|
| 11937994 | 11.9 x 16.8 x 17.5 | 1 |

Parafilm™ is a registered trademark of Amcor PLC. Tough-Tags™ is a registered trademark of Diversified Biotech, Inc.




Fisherbrand™ Label Station

Organise your labels for easy accessibility

Portable label station holds up to six standard size boxes of Tough-Tags™, Cryo-Babies™, and Tough-Spots™

Just drop the box in the label station and feed the paper through the slot. Labels shown not included, order separately

- Acrylic station includes a writing surface, a cutting edge for easy dispensing and non-skid rubber feet
- Counter-sunk holes provided for wall mounting
- Translucent blue

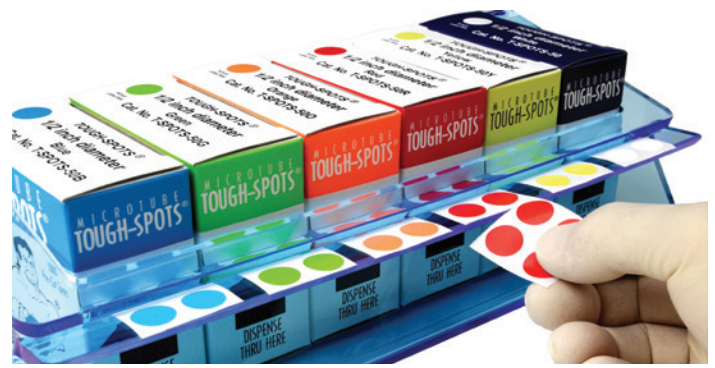

| Cat. No. | Colour | L x W x H mm | Qty |
|----------|--------|---------------|-----|
| 11947994 | Blue | 28 x 12 x 8.9 | 1 |

Tough-Tags™, Cryo-Babies™ and Tough-Spots™ is a registered trademark of Diversified Biotech, Inc.



Unique safety lid

Fisherbrand™ Mini Bin Foil Dispenser

Safely cut foil to any size

The mini bin foil dispenser is a great alternative to wasteful giant-sized foil rolls and expensive pre-cut foil squares

- Compact to fit in limited work space areas and easy to transport
- Safety guard feature to prevent injuries
- Foam strip prevents foil retracting back into unit
- For use with 102 mm or 153 mm length foil rolls (sold separately)
- Made from ABS

| Cat. No. | L x W x H cm | Qty |
|----------|----------------|-----|
| 11957994 | 19 x 10 x 11.2 | 1 |


Fisherbrand™ Foil Rolls

Designed for easy use and minimal waste

- Custom-sized pieces of foil can easily be created for covering small surfaces such as beakers, flasks, tubes and Petri dishes
- Ideal for autoclaving, freezing, incubating, protecting labware, lining, sealing and general storage
- Two sizes available. designed to fit the mini bin foil dispenser



| Cat. No. | Roll Width cm | Roll Length m | Qty |
|----------|---------------|---------------|-----|
| 11967994 | 10.2 | 152.4 | 1 |
| 11977994 | 15.3 | 152.4 | |

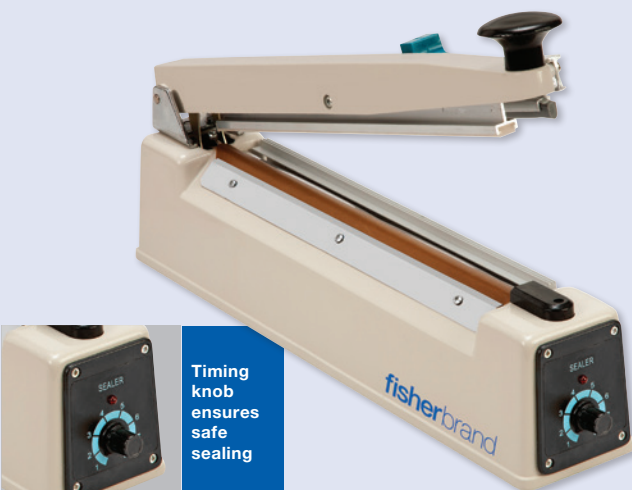

Timing knob ensures safe sealing

Fisherbrand™ Heat Sealer with Cutter

Quickly and efficiently seal polyethylene or polypropylene bags, sheets or flat tubing from 0.01 to 0.8 mm thick

Sliding blade removes excess material 19 mm from the seal

- Can seal a bag up to 12 in (30.5 cm) long
- Heating element measures 31 cm long and is only hot when sealing
- Includes one replacement part kit
- Additional replacement part kits can be ordered separately

| Cat. No. | Power Supply | D x W x H cm | Qty |
|----------|-----------------------------------|-------------------|-----|
| 15350360 | 220/230 VAC, U.K | 8.6 x 44.6 x 16.6 | 1 |
| 15360360 | 220/230 VAC, Cont. Europe/Israel | | 1 |
| 15350407 | Heat Sealer Replacement Parts Kit | | 1 |

Replacement Parts Kit contains: Blades (2), Heating element (2), PTFE cloths (2)



Distributed by Fisher Scientific. Contact us today:

Austria: fishersci.at **Belgium:** fishersci.be **Denmark:** fishersci.dk
Germany: fishersci.de **Ireland:** fishersci.ie **Italy:** fishersci.it
Finland: fishersci.fi **France:** fishersci.fr **Netherlands:** fishersci.nl
Norway: fishersci.no **Portugal:** fishersci.pt **Spain:** fishersci.es
Sweden: fishersci.se **Switzerland:** fishersci.ch **UK:** fishersci.co.uk

© 2023 Thermo Fisher Scientific Inc. All rights reserved.
Trademarks used are owned as indicated at fishersci.com/trademarks.

 **fisher scientific**
part of Thermo Fisher Scientific